Oral-B®

Professional Care Center

OC 15 525 A

powered
by **BRAUN**

Type 3718


Deutsch	4, 9, 94
English	10, 15, 94
Français	16, 21, 94
Español	22, 27, 94
Português	28, 33, 94
Italiano	34, 39, 94
Nederlands	40, 45, 94
Dansk	46, 51, 94
Norsk	52, 57, 94
Svenska	58, 63, 94
Suomi	64, 69, 94
Ελληνικά	70, 75, 94
Türkçe	76, 94
Русский	81, 86, 94
عربي	93, 87, 94

Internet:

www.oralb.com
www.braun.com

- (D)** Haben Sie Fragen zu diesem Produkt?
- (A)** Rufen Sie an (gebührenfrei):
00 800 27 28 64 63
- (CH)** **08 44 - 88 40 10**
- (GB)** If you should require further assistance, please call our Consumer Service line.
0 800 73 11 792 (free call)
- (IRL)** **Helpline**
1 800 509 448 (free call)
- (F)** Appelez le Service Consommateurs Oral-B :
0 810 365 855 (numéro azur)
- (B)** **02-711 92 11**
- (LU)** Vous avez des questions sur ce produit? Appelez:
4 00 50 51
- (E)** Servicio al consumidor para España:
901 11 61 84
- (P)** Serviço ao Consumidor para Portugal:
808 20 00 33
- (I)** Servizio consumatori:
02-667 86 23
- (NL)** Heeft u vragen over dit product? Bel Braun/Oral-B Consumenten-ijfolijn:
0800-4 45 53 88
- (DK)** Har du spørgsmål om produktet? Så ring
70 15 00 13
- (N)** Spørsmål om dette produktet? Ring
22-72-88-10
- (S)** Frågor om apparaten? Ring Kundservice
020 - 21 33 21
- (FIN)** Onko Sinulla kysyttävää tuotteesta? Soita
09 452 872
- (GR)** **Τηλέφωνο για πληροφορίες**
210 9478700
- (TR)** **Danışma Hattı**
0212-473-75-85
- (RUS)** **Телефонный помощь**
8-800-200-1444

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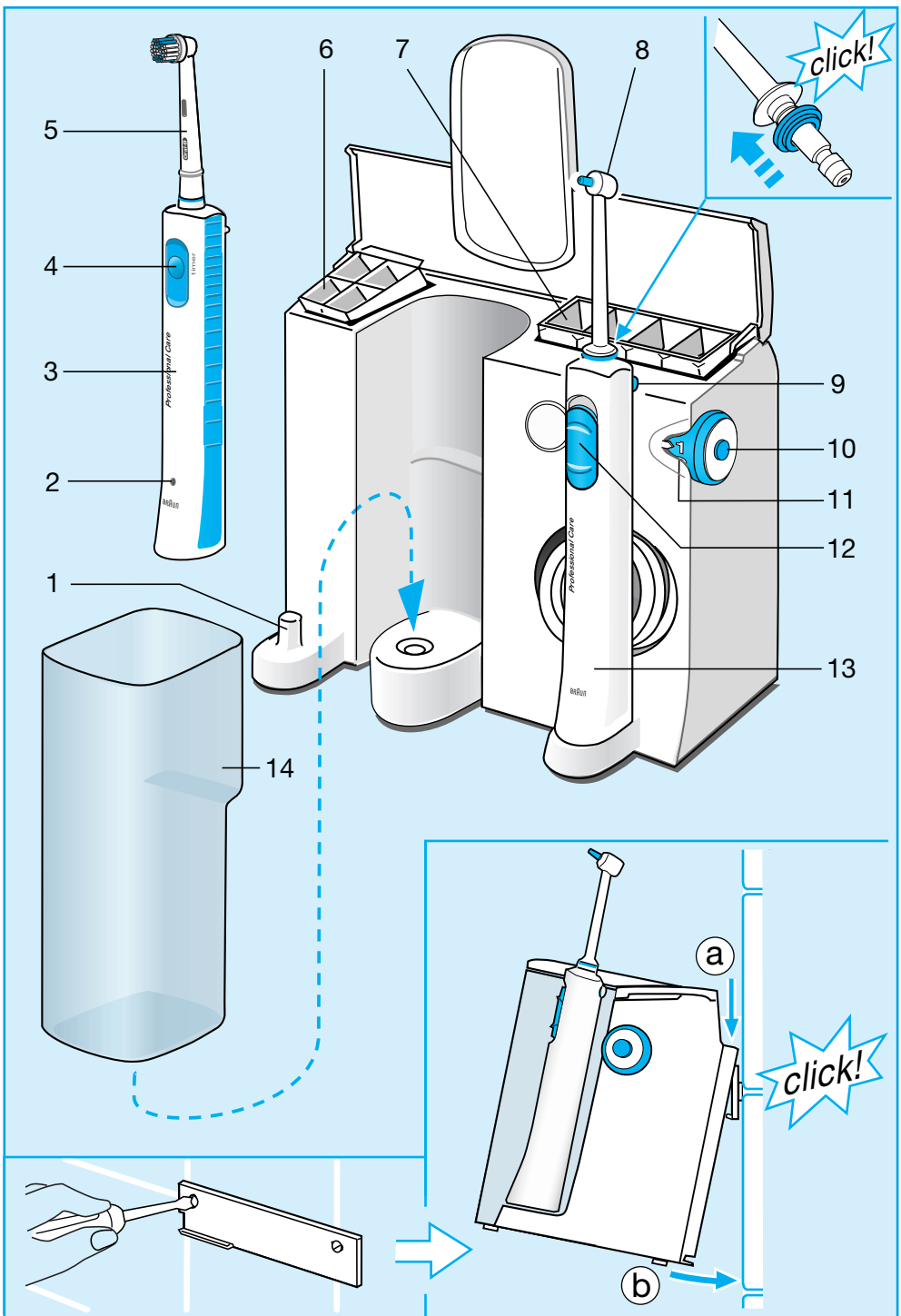The Oral-B Professional Care Center has been developed together with leading dental professionals to effectively remove plaque from all areas of your teeth and to attack plaque bacteria.

The Oral-B Professional Care toothbrush adds a gentle pulsating movement to its oscillation to create a deep-cleaning 3D brushing action. The toothbrush should be used twice a day for at least two minutes, the minimum brushing time recommended by dental professionals.

The irrigator removes remaining food particles from between your teeth, while it also effectively massages your gums. Always use the water jet after brushing your teeth

Please read the use instructions carefully before first use.

Important

The appliance must be mounted in a way that it cannot fall into water.

The plug should always be accessible. Periodically check the cord for damage. A damaged or non-functioning appliance should no longer be used. Take it to a Braun Service Centre.

Description

Toothbrush

- ① Charging base
- ② Charge indicator
- ③ Toothbrush handle (rechargeable)
- ④ Switch
- ⑤ Brushhead
- ⑥ Brushhead compartment

Irrigator

- ⑦ Jet compartment
- ⑧ Irrigator jet
- ⑨ Jet release button
- ⑩ On/off switch
- ⑪ Water pressure control dial
- ⑫ Sliding switch (water stop and go)
- ⑬ Irrigator handle
- ⑭ Water container (600 ml/20 fl. oz.)

Specifications

Power supply: see bottom of housing

Power intake: 24 W

Using the toothbrush

Connecting and charging

Plug the appliance into an electrical outlet. Place the handle ③ on the charging base ①, the charge indicator ② will light up.

A full charge will take at least 16 hours and will provide 30 minutes operation time. For everyday use, the handle can be stored on the charging base to maintain it at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, the following procedure should be done every 6 months: Do not place the handle back onto the charging base in order to discharge it by regular use. Once discharged, recharge the handle on the charging base.







Operating the switch

Press the switch ④ to operate the toothbrush, press again to stop it.

Removing plaque

- The toothbrush can be used with any leading brand of toothpaste.
- To avoid splashing, guide the brushhead to your teeth before switching on.
- Guide the brushhead slowly from tooth to tooth. Hold the brushhead in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces. Do not press too hard nor scrub, simply let the brush do all the work.

During the first few days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or hygienist.

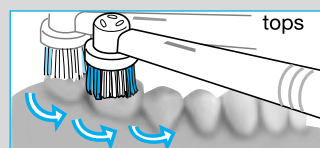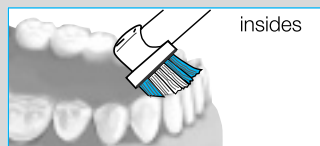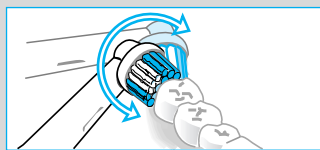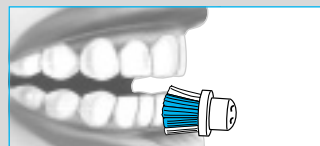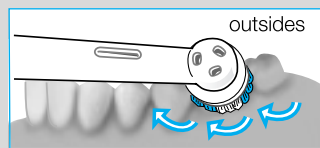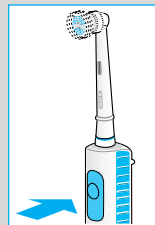
Pressure sensor

For safety and optimum brushing results, the toothbrush has a pressure control system: If too much pressure is applied, the oscillating movement of the brushhead will continue, but the pulsation will stop. There is a noticeable change in feeling and sound when the pressure sensor is activated.

Memory timer

An incorporated timer memorises the elapsed brushing time, even when the handle is briefly switched off (max. 30 seconds). After two minutes – which is the minimum recommended brushing time – the timer signals by a short stutter.

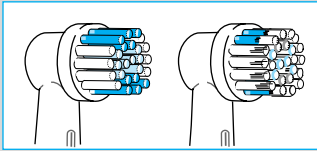
INDICATOR® bristles

The INDICATOR® bristles help you monitor and improve the brushing performance.

With proper brushing, twice a day for two minutes, the blue colour on the INDICATOR® bristles will disappear half-way down within 3-4 months, signalling it is time to replace the brushhead.

Replacement brushheads

Available at your retail store or Braun Service Centres.

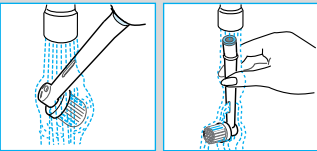
After use

Rinse the brushhead thoroughly under running water for several seconds with the handle switched on. Switch off the handle and take off the brushhead. Clean both parts separately under running water, then wipe them dry.

Environmental notice

The toothbrush handle is provided with a nickel-hydride rechargeable battery. To protect the environment, do not dispose of the product in the household waste at the end of its useful life. Open the handle as shown, remove the battery and dispose of it at your retail store or at appropriate collection sites according to national or local regulations.

Caution: Opening the handle will definitely destroy it.

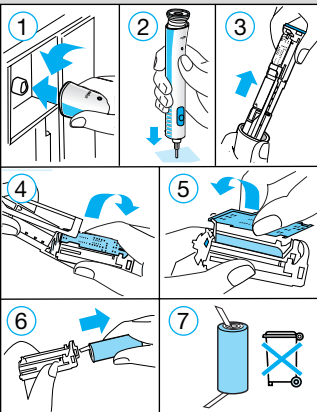Using the Irrigator

Before using for the first time, fill the water container and allow the water to run through completely by holding the irrigator handle over a washbasin.

1. Attach one of the irrigator jets ⑧ to the handle ⑬ until it locks in place. Fill the water container with water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch ⑩. (When doing so, the sliding switch ⑫ on the handle should not be pushed up.)
3. Select the water pressure with the pressure control dial ⑪ (1 = gentle, 5 = strong). We recommend that you start with setting (1), especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over the washbasin and guide the jet to your teeth and gums. Slide the handle switch up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.

The irrigator jet offers two different functions depending on the position of the jet tip:

Single jet:





Jet tip in position (a) to clean the teeth and the interdental spaces and remove resistant food particles

Multi jet:


Jet tip in position (b) to massage the gums and keep them healthy.

First, turn off the water stream by sliding the handle switch down (to avoid unintended water spraying). The jet can easily be switched from one function to the other by pressing the jet tip in or out. Begin with the single jet. Guide the water stream to the teeth and especially to the interdental spaces. Then massage the gums with the multi jet.

To ensure the jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

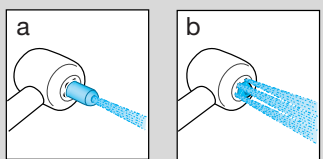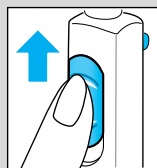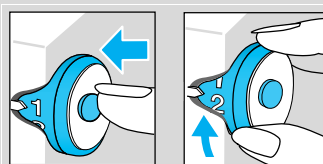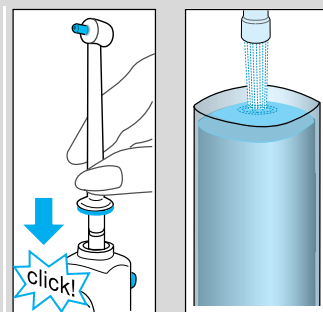Max. operating time: 15 minutes
Cooling off time: 2 hours

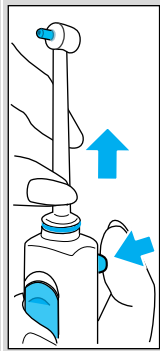After use

Always empty the water container entirely. Otherwise bacteria could collect in the remaining water, if the unit is not used for some time. Should your dentist have recommended use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use with chlorhexidine, to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator.

To take off the irrigator jet, press the jet release button. The jets can be stored in the jet compartment. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth. The water container, the brushhead compartment and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

At the end of the product's useful life, please dispose of it at appropriate collection points provided in your country.

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

Garantie/Kundendienststellen**Guarantee and Service Centers****Bureaux de garantie et centrales service après-vente****Oficinas de garantía y oficinas centrales del servicio****Entidade de garantia e centros de serviço****Uffici e sede centrale del servizio clienti****Garantiebureaux en service-centrales****Garantikontorer og centrale serviceafdelinger****Garanti og servicecenter****Servicekontor och centrala kundtjänstplatser****Asiakaspalvelu ja tekninen neuvonta****Εγγύηση και Εργαστήριο επισκευών****Garanti bürolari ve merkezi servis yerleri****Гарантийные и сервисные центры****مراكز الكفالات والصيانة****Deutschland**

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Hong Kong

Audio Supplies Company,
Room 506,
St. George Building,
2 ICE House Street,
Hong Kong,
☎ (852) 2524 9377

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Tehran Bouran Company
145 Ghaem Magham Farahani
Tehran
☎ +982 18 31 27 66

Ireland (Republic of)

Gavin's Electronics,
83-84, Lower Camden Street,
Dublin 2,
☎ 1800 509 448

Israel

S. Schestowitz Ltd.,
8 Shacham Str.,
Tel-Aviv, 49517,
☎ 1 800 335 959

Italia

Servizio Consumatori Braun
Gillette Group Italy S.p.A.,
Via G.B. Pirelli, 18,
20124 Milano,
☎ 02 / 6678623

Jordan

Interbrands
Wadi Al-Sir, Adwiat alhekma
St. Mekadabi bld.
Amman
☎ +692 582 75 67

Kenya

Radbone-Clark Kenya Ltd.,
P.O. Box 40833,
Mombasa Road, Nairobi,
☎ 2 82 12 76

Korea

Gillette Korea Ltd.
144-27 Samsung-dong,
Kangnam-ku, Seoul, Korea,
☎ 080-920-6000

Kuwait

Union Trading Company,
Braun Service Center,
P.O. Box 28 Safat,
Safat Code 13001, Kuwait,
☎ 04 83 32 74