
BRAUN

Oral-B

MD 15

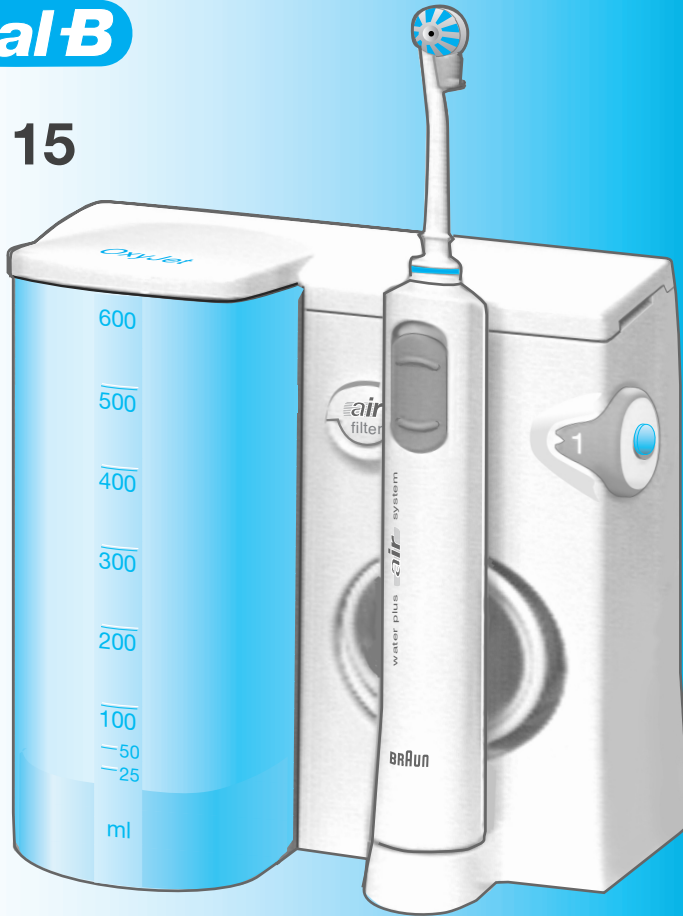OxyJet®

Type 4715

Braun Infoline

Deutsch	4, 36, 45
English	6, 36, 45
Français	8, 36, 45
Español	10, 36, 45
Português	12, 37, 45
Italiano	14, 37, 45
Nederlands	16, 37, 45
Dansk	18, 38, 45
Norsk	20, 38, 45
Svenska	22, 38, 45
Suomi	24, 38, 45
Ελληνικά	26, 39, 45
Türkçe	29, 45
Русский	31, 39, 45
عربي	36, 40, 45

- (D)** Haben Sie Fragen zu diesem Produkt?
Rufen Sie an: (in Deutschland und Österreich zum Nulltarif)
00 800 27 28 64 63
00 800 BRAUNINFOLINE
- (CH)** **08 44 - 88 40 10**
- (E)** Servicio al consumidor para España:
901-11 61 84
- (P)** Serviço ao Consumidor para Portugal:
080 820 00 33
- (I)** Servizio consumatori:
(02) 6 67 86 23
- (NL)** Heeft u vragen over dit produkt?
Bel Braun Consumenten-infolijn:
(070) 4 13 16 58
- (B)** Vous avez des questions sur ce produit ?
Appelez Braun Belgique
(02) 711 92 11
- (DK)** Har du spørgsmål om produktet? Så ring
70 15 00 13
- (N)** Spørsmål om dette produktet? Ring
88 02 55 03
- (S)** Frågor om apparaten? Ring Kundservice
020 - 21 33 21
- (FIN)** Onko Sinulla kysyttävää tuotteesta? Soita
0203 77877
- (TR)** Türkiye'deki Tüketici Danışma Servisi:
0800 261 19 53

Internet:
www.braun.com
www.oralb.com

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English

The Braun Oral-B OxyJet® irrigator mixes air into water in the form of micro-bubbles. These micro-bubbles in the water stream attack plaque bacteria.

At the same time, it removes remaining food particles from between your teeth and massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

The plug should always be accessible. Periodically check the cord for damage. Should this be the case, take the appliance to a Braun Service Centre. A damaged or non-functioning unit should no longer be used.

Description

- a Irrigator jet
- b Jet switch
- c Jet compartment
- d Jet release button
- e On/off switch
- f Water pressure control dial
- g Sliding switch (water stop and go)
- h Irrigator handle
- i Water container (600 ml/20 fl. oz.)

Specifications

Power supply: see bottom of housing

Power intake: 24 W

Using the OxyJet® Irrigator


Before using for the first time, fill the water container and allow the water to run through completely by holding the irrigator handle over a washbasin.

1. Attach one of the irrigator jets to the handle until it locks in place. Fill the water container with water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch. (When doing so, the sliding switch on the handle should not be pushed up.)
3. Select the water pressure with the pressure control dial (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over the wash-basin and guide the jet to your teeth and gums. Slide the handle switch up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.


You can choose between two different jet functions by sliding the jet switch down or up (before doing so, turn off the water stream by pushing the sliding switch down):

Rotating stream: for general cleaning and gum massage (jet switch in lower position)

The rotation of the water stream and its micro-bubbles is designed for easy penetration below the gumline.

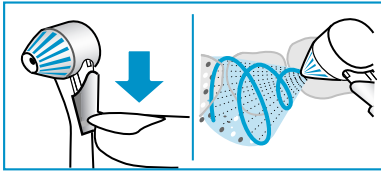Mono stream: straight water stream for site-specific cleaning (jet switch in upper position)

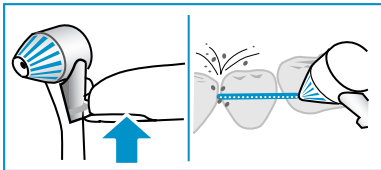Max. operating time: 15 minutes
Cooling off time: 2 hours

After use

Always empty the water container entirely. Otherwise bacteria could collect in the remaining water, if the unit is not used for some time. If your dentist has recommended use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator. The irrigator will automatically turn off after 10 minutes if you forget to turn it off manually. In this case, to turn it on again, you have to press the on/off switch twice.

To take off the irrigator jet, press the jet release button. The jets can be stored in the jet compartment. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth.

The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the EMC-Requirements as laid down by the Council Directive 89/336/EEC and to the Low Voltage Regulation (73/23 EEC).

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English

Guarantee

We grant 1 year guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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Danmark A/S,