**Oral-B®**

*ProfessionalCare*  
*6500 WaterJet*

MD 16

powered  
by **BRAUN**

Deutsch	4, 21
English	6, 21
Français	8, 21
Español	10, 22
Português	12, 22
Italiano	14, 22
Nederlands	16, 23
Ελληνικά	18, 23

**DE AT 00 800 27 28 64 63**

**CH 08 44 - 88 40 10**

**GB 0 800 731 1792**

**IE 1 800 509 448**

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# English

The **Oral-B ProfessionalCare WaterJet** removes remaining food particles from between your teeth, while it also effectively massages your gums. Always use the WaterJet after brushing your teeth.

Please read the use instructions carefully before first use.

## Important

- The appliance must be mounted in a way that it cannot fall into water or drop to the floor. The plug should always be accessible.
- Periodically check the cord for damage. A damaged or non-functioning unit should no longer be used. If the cord is damaged, take the charging unit to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age 12.
- This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are supervised by a person responsible for their safety.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not place the appliance in water or other liquid.
- Do not reach for an appliance that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Braun Service Centre. Do not disassemble the product.
- Do not insert any object into any opening of the appliance.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments or charges not recommended by the manufacturer.

- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

## Description

- Ⓐ Irrigator jet
- Ⓑ Jet compartment
- Ⓒ Jet release button
- Ⓓ On/off switch
- Ⓔ Water pressure control dial
- Ⓕ Handle switch (water stop and go)
- Ⓖ Irrigator handle
- Ⓗ Water container (600 ml)
- Ⓘ Cord storage

## Specifications:

For voltage specifications please refer to the bottom of the charging unit.

## Using the irrigator


Before using for the first time, fill the water container **Ⓗ** and allow the water to run through completely by holding the irrigator handle over a washbasin.

1. Attach one of the irrigator jets **Ⓐ** to the handle **Ⓖ** until it locks into place. Fill the water container with luke warm water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch **Ⓓ**. (When doing so, the handle switch **Ⓕ** should be in the lower position = water stop).
3. Select the water pressure with the pressure control dial **Ⓔ** (1 = gentle, 5 = strong). We recommend that you start with setting (1), especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.

The irrigator jet offers two different functions depending on the position of the jet tip:

Single jet:

Jet tip in position (a) to clean the teeth and the interdental spaces and remove resistant food particles.

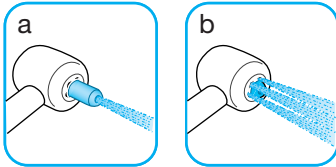Multi jet:


Jet tip in position (b) to massage the gums and keep them healthy.

First, turn off the water stream by sliding the handle switch down (to avoid unintended water spraying). The jet can easily be switched from one function to the other by pressing the jet tip in or out. Begin with the single jet. Guide the water stream to the teeth and especially to the interdental spaces. Then massage the gums with the multi jet.

Max. operating time: 15 minutes

Cooling off time: 2 hours

To ensure the irrigator jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

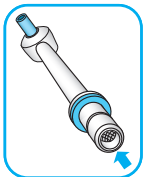**After use**

Always empty the water container entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended the use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator. To take off the irrigator jet, press the jet release button (ⓐ). The jets can be stored in the jet compartment (ⓑ). Wipe the handle dry and place it back onto its base.

**Cleaning**

Regularly clean the housing with a damp cloth. The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives EMC 2004/108/EC and Low Voltage 2006/95/EC.

Please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

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## English

### Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country

where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Customer Service Centre.

### For UK only:

This guarantee in no way affects your rights under statutory law.

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Munddusche mit dem Ein-/Aus-Schalter ausschalten. Entriegelungsknopf © drücken, um die Aufsteckdüse abzunehmen. Wischen Sie das Handstück ab und setzen Sie es zurück auf seinen Sockel. Die Aufsteckdüsen können im Aufbewahrungsfach ⑤ aufbewahrt werden.