Oral-B®

ProfessionalCare

8000 OxyJet

powered by **BRAUN**

Deutsch	4	D	A	00 800 27 28 64 63
English	7	CH		08 44 - 88 40 10
Français	10	UK		0 800 416 434
Español	13	IRL		1 800 509 448
Português	16	F		0 810 365 855
Italiano	19	BE		0 800 14 592
Nederlands	22	L		4 00 50 51
Ελληνικά	25	ES		901 11 61 84
		PT		808 20 00 33
		I		800 124 600
		NL		0800-4 45 53 88
		GR		800 801 3456

Internet:

www.oralb.com

www.braun.com

www.service.braun.com

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English

The **Oral-B ProfessionalCare 8000 OxyJet** irrigator mixes air into water in the form of micro-bubbles. It helps refresh your mouth. At the same time, the irrigator removes remaining food particles from between your teeth and massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

- The appliance must be mounted in a way that it cannot fall into water. The plug should always be accessible.
- Periodically check the cord for damage. A damaged or non-functioning appliance shall no longer be used. Take it to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are given supervision by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.

Description

- Ⓐ Irrigator jet
- Ⓑ Jet switch
- Ⓒ Jet compartment
- Ⓓ Jet release button
- Ⓔ On/off switch
- Ⓕ Water pressure control dial
- Ⓖ Handle switch (water stop and go)
- Ⓗ Irrigator handle

- Ⓘ Water container (600 ml)
- Ⓚ Cord storage

Specifications

Power supply: see bottom of housing

Power intake: 24 W


Using the irrigator

Before using for the first time, fill the water container ① and allow the water to run through completely by holding the irrigator handle ② over a washbasin.


1. Attach one of the irrigator jets ③ to the handle until it locks into place. Fill the water container with water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch ④.
(When doing so, the handle switch should be in the lower position = water stop).
3. Select the water pressure with the pressure control dial ⑤ (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch ⑥ up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.

You can choose between two different jet functions by sliding the jet switch ⑦ down or up (before doing so, turn off the water stream by sliding the handle switch down):

Rotating stream: for general cleaning and gum massage. The rotation of the water stream and its micro-bubbles is designed for easy penetration (jet switch in lower position).

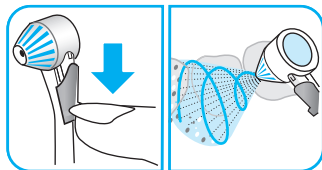
Mono stream: straight water stream for site-specific cleaning (jet switch in upper position).

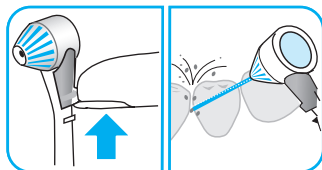Max. operating time: 15 minutes

Cooling off time: 2 hours

To ensure the irrigator jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

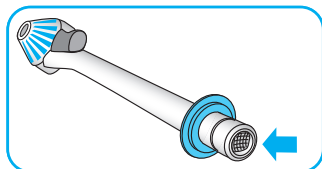After use

Always empty the water container entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended the use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator. The irrigator will automatically turn off after 10 minutes if you forget to turn it off manually. In this case, to turn it on again, you have to press the on/off switch twice.

To take off the irrigator jet, press the jet release button Ⓐ. The jets can be stored in the jet compartment Ⓒ. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth. The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives EMC 2004/108/EC and Low Voltage 2006/95/EC.

Please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brushheads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

For UK only:

This guarantee in no way affects your statutory rights.