Oral-B®

*Professional
Care OxyJet®*

MD 17

powered
by **BRAUN**

Type 3719

Deutsch	4, 44, 51
English	6, 44, 51
Français	8, 44, 51
Español	10, 44, 51
Português	12, 45, 51
Italiano	14, 45, 51
Nederlands	16, 45, 51
Dansk	18, 46, 51
Norsk	20, 46, 51
Svenska	22, 46, 51
Suomi	24, 46, 51
Polski	26, 47, 51
Český	28, 47, 51
Slovenský	30, 48, 51
Magyar	32, 48, 51
Ελληνικά	34, 48, 51
Türkçe	37, 51
Русский	39, 49, 51
عربي	43, 50, 51

Internet:

www.oralb.com
www.braun.com

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English

The Oral-B Professional Care OxyJet® irrigator mixes air into water in the form of micro-bubbles. These micro-bubbles in the water stream attack plaque bacteria. At the same time, it removes remaining food particles from between your teeth and massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

The appliance must be mounted in a way that it cannot fall into water.

The plug should always be accessible.

Periodically check the cord for damage.

A damaged or non-functioning unit should no longer be used. Take it to a Braun Service Centre.

Description

- Ⓐ Irrigator jet
- Ⓑ Jet switch
- Ⓒ Jet compartment
- Ⓓ Jet release button
- Ⓔ On/off switch
- Ⓕ Water pressure control dial
- Ⓖ Handle switch (water stop and go)
- Ⓗ Irrigator handle
- Ⓘ Water container (600 ml/20 fl. oz.)

Specifications

Power supply: see bottom of housing

Power intake: 24 W

Using the OxyJet® Irrigator


Before using for the first time, fill the water container and allow the water to run through completely by holding the irrigator handle over a washbasin.

1. Attach one of the irrigator jets to the handle until it locks in place. Fill the water container with water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch. (When doing so, the handle switch should be in the lower position = water stop.)
3. Select the water pressure with the pressure control dial (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over the washbasin and guide the jet to your teeth and gums. Slide the handle switch up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.


You can choose between two different jet functions by sliding the jet switch down or up (before doing so, turn off the water stream by pushing the handle switch down):

Rotating stream: for general cleaning and gum massage (jet switch in lower position)

The rotation of the water stream and its micro-bubbles is designed for easy penetration below the gumline.

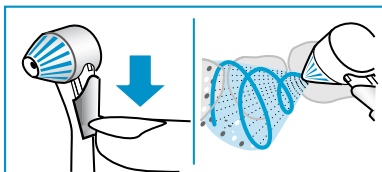Mono stream: straight water stream for site-specific cleaning (jet switch in upper position)

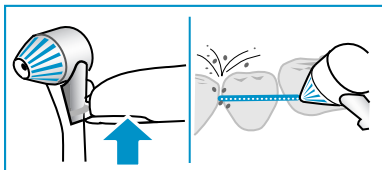Max. operating time: 15 minutes
Cooling off time: 2 hours

To ensure the irrigator jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

After use

Always empty the water container entirely to prevent bacteria from collecting in the stand-ing water. If your dentist has recommended use of a chlorhexidine solution with your ir-rigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator. The irrigator will automatically turn off after 10 minutes if you forget to turn it off manually. In this case, to turn it on again, you have to press the on/off switch twice.

To take off the irrigator jet, press the jet release button. The jets can be stored in the jet compartment. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth.

The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

At the end of the product's useful life, please dispose of it at appropriate collection points provided in your country.

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English

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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Gillette (Shanghai) Sales Co. Ltd.
Gillette de Colombia S.A.,
Gillette Group Danmark A/S,