
PROFESSIONAL
CARE®

powered by **BRAUN**

WATERJET

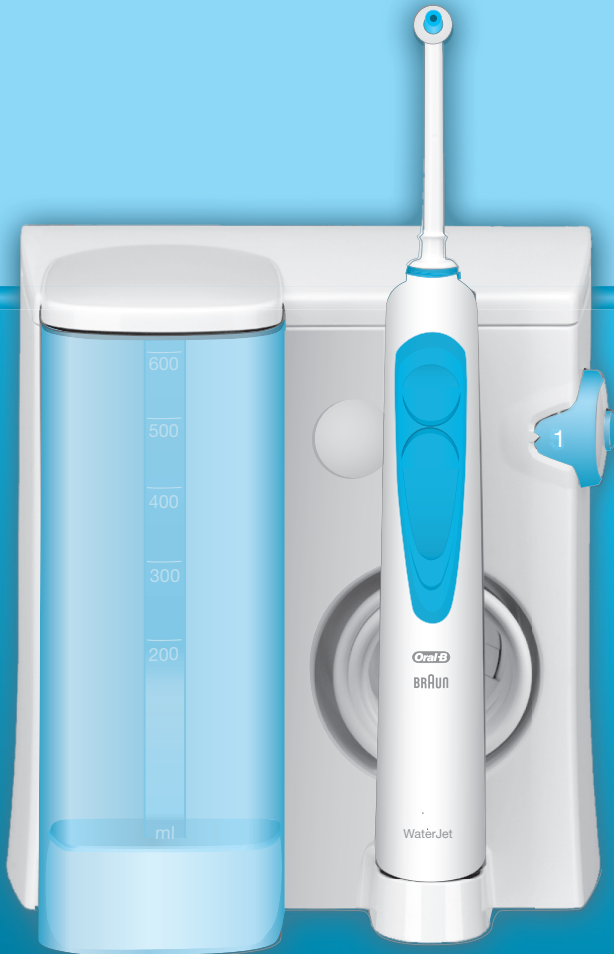Deutsch	4	DE AT	00 800 27 28 64 63
English	7	CH	08 44 - 88 40 10
Français	9	GB	0 800 731 1792
Español	12	IE	1 800 509 448
Português	15	FR	0 800 944 803 (appel gratuit depuis un poste fixe)
Italiano	18	BE	0 800 14 592
Nederlands	20	LU	800 21172
Ελληνικά	23	ES	901 11 61 84
		PT	808 20 00 33
		IT	800 124 600
		NL	0800-4 45 53 88
		GR	800 801 3457

Internet:

www.oralb.com

www.braun.com

www.service.braun.com

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English

The Oral-B ProfessionalCare Waterjet provides a fast, yet gentle removal of food particles between your teeth while also stimulating and massages your gums. Used in combination with proper tooth brushing, it helps improve your gum health. For most effective results always use the Waterjet after tooth brushing.

Please read the usage instructions carefully before first use.

IMPORTANT

- The appliance must be mounted in a way that it cannot fall into water or drop to the floor. The plug should always be accessible.
- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- This product is not intended for use by children under age 12. Persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can use the product, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Oral irrigators shall be used only with liquids defined in the manufacturer's instructions.
- Use this product only for its intended use as described in this manual. Do not use attachments or chargers which are not recommended by the manufacturer.

WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.

- Do not place the appliance in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- Do not disassemble the product.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

- a Irrigator jet
- b Jet compartment
- c Jet release button
- d On/off switch
- e Water pressure control dial
- f Handle switch (water stop and go)
- g Irrigator handle
- h Water container (600 ml)
- i Cord storage

Specifications

For voltage specifications please refer to the bottom of the charging unit.

Using the irrigator

Before using for the first time, fill the water container (h) and allow the water to run through completely by holding the irrigator handle (g) over a washbasin.

1. Attach one of the irrigator jets (a) to the handle until it locks into place (1). Fill the water container with lukewarm water. If you want to add mouthwash into the water, please make sure you fill the water container with lukewarm water first, then add the mouthwash.
2. Turn on the irrigator by pressing the on/off switch (d) (2). When doing so, the handle switch (f) should be in the lower position = water stop).
3. Select the water pressure with the pressure control dial (e) (1 = gentle, 5 = strong) (3). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch up to start the water stream (4). Keep your mouth slightly open to allow the water to run into the washbasin.

The irrigator jet offers two different functions depending on the position of the jet tip:

Single jet:

Jet tip is in position (6) to clean teeth and the interdental spaces and remove resistant food particles.

Multi jet:

Jet tip is in position (7) to massage the gums and keep them healthy.

First, turn on the water stream by sliding the handle switch down (to avoid unintended water spraying). The jet can easily be switched from one function to the other by pressing the jet tip in or out. Begin with the single jet. Guide the water stream to the teeth; especially to the interdental spaces. Then massage the gums with the multi jet.

Maximal operation time: 15 minutes

Cooling off time: 2 hours

In general we recommend to replace the irrigator jet every 3 months.

After use

Always empty the water container entirely. If your dentist has recommended the use of mouthwash with your irrigator, allow some clean water to run through the irrigator before you empty the container entirely to prevent any deposit build-up. Make sure you keep the irrigator switched on until the water stops running through the irrigator jet, to prevent bacteria from collecting in the standing water.

Press the on/off switch to turn off the irrigator.

To take off the irrigator jet, press the jet release button (c) (5). Always wipe the handle and the irrigator jet dry; place the handle back onto its base and store the irrigator jet in the jet compartment (b).

Cleaning

Regularly clean the housing with a damp cloth. The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

Please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

For UK only:

This guarantee in no way affects your statutory rights.

30-Days Money Back Guarantee

Try Oral-B ProfessionalCare Waterjet irrigator starting from the day of purchase. If you are not satisfied with this product, please return your handle, refills and charging unit in their original box with your original cash register receipt to the address below within 30 days of purchase.

Braun Cashback
PO Box 504
Leicester
LE940AE

Terms and conditions apply, for further details, please contact our customer service.

Warranty replacement brush heads and accessories

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads and/or accessories.

Oral-B does not recommend the use of non Oral-B replacement brush heads and/or accessories.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads and/or accessories. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads and/or accessories, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads and/or accessories.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads and/or accessories on handle wear.

All Oral-B replacement brush heads and/or accessories carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads and/or accessories or handle parts under any other brand name.