Oral-B®

Professional Care

6500 WaterJet Center

OC 16.525

powered by **BRAUN**

Deutsch	4	DE AT	00 800 27 28 64 63
English	10	CH	08 44 - 88 40 10
Français	16	GB	0 800 731 1792
Español	22	IE	1 800 509 448
Português	28	FR	0 800 944 803 (appel gratuit depuis un poste fixe)
Italiano	34	BE	0 800 14 592
Nederlands	40	LU	4 00 50 51
Ελληνικά	46	ES	901 11 61 84
		PT	808 20 00 33
		IT	800 124 600
		NL	0800-4 45 53 88
		GR	800 801 3457

Internet:

www.oralb.com
www.braun.com
www.service.braun.com

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The **Oral-B ProfessionalCare Center** has been developed together with leading dental professionals to effectively remove plaque from all areas of your teeth and to attack plaque bacteria.

The **Oral-B ProfessionalCare** toothbrush provides outstanding cleaning for clean teeth and healthy gums.

The **Oral-B ProfessionalCare** irrigator removes remaining food particles from between your teeth, while it also effectively massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

- The appliance must be mounted in a way that it cannot fall into water or drop to the floor. The plug should always be accessible.
- Periodically check the cord for damage. A damaged or non-functioning unit should no longer be used. If the cord is damaged, take the charging unit to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age 12.
- This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by a person responsible for their safety.
- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not place the appliance in water or other liquid.
- Do not reach for an appliance that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Braun Service Centre. Do not disassemble the product except when disposing of the battery.
- When taking out the battery from the toothbrush handle for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Do not insert any object into any opening of the appliance.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments or charges not recommended by the manufacturer.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

Toothbrush

- ① Charger base
- ② Charge indicator
- ③ Toothbrush handle (rechargeable)
- ④ On/off button
- ⑤ Brush head
- ⑥ Brush head compartment

Irrigator

- ⑦ Jet compartment
- ⑧ Irrigator jet
- ⑨ Jet release button
- ⑩ On/off switch
- ⑪ Water pressure control dial
- ⑫ Handle switch (water stop and go)
- ⑬ Irrigator handle
- ⑭ Water container (600 ml)
- ⑮ Cord storage

Specifications:


For voltage specifications please refer to the bottom of the charging unit.





Using your toothbrush

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and is designed for use in the bathroom.

- Plug the appliance into an electrical outlet and place the toothbrush handle ③ on the charger base ①.
- The **green** charge indicator ② light flashes while the handle is being charged.
- Once the handle is fully charged, the green charge indicator light shines continuously for 5 seconds and then it turns off. A full charge takes at least 18 hours and enables approximately one week of regular brushing (twice a day, 2 minutes). To check the charging status when the charging light is off, shortly lift up and return your handle on the charging unit. The charge indicator light will shine permanently for 5 seconds and turn off when your handle is fully charged.
- The **red** charge indicator ② light starts blinking when the rechargeable battery is low. Now it is time to recharge your toothbrush.
- For everyday use, the toothbrush handle can be stored on the charger base to maintain it at full power. Overcharging is impossible. However, for environmental reasons, Braun Oral-B recommends to unplug the charging unit until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, refrain from placing the handle on the charger base and fully discharge it through regular use at least every 6 months.

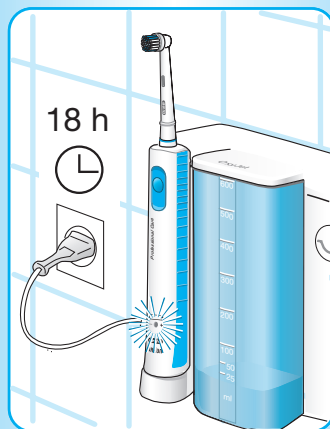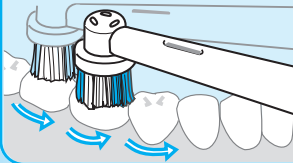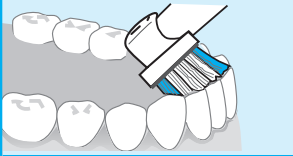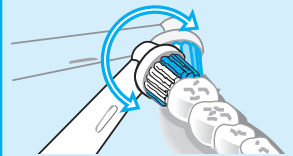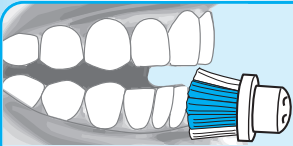
Brushing technique

- Attach one of the brush heads to the handle. They can be used with any kind of toothpaste.
- To avoid splashing, guide the brush head to your teeth before switching the appliance on.
- Guide the brush head slowly from tooth to tooth.
- Hold the brush head in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.
- To thoroughly remove plaque, brush for at least two minutes.

Pressure Sensor

For optimal brushing results **Oral-B ProfessionalCare** has a pressure control feature to help you brush gently. This helps protect your gums from excessive brushing and helps make your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition to feeling it, you can also hear a different sound.

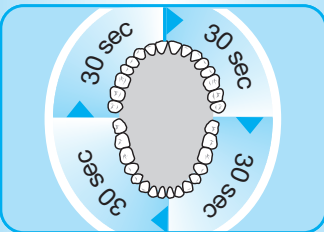Professional Timer


A short stuttering sound at 30-second intervals, which you can hear and feel, reminds you to brush the four quadrants of your mouth equally for a thorough consistent clean. A long stuttering sound indicates when the dentist recommended 2-minute brushing time has been reached.

The timer memorises the elapsed brushing time, even when the handle is briefly switched off during brushing. Only when pausing longer than 30 sec. the timer resets to the initial starting point.

In the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or dental hygienist.

Dentists recommend that you replace any toothbrush every 3 months. The INDICATOR bristles help you monitor and improve your brushing performance. With proper brushing, twice a day for two minutes, the blue colour on the INDICATOR bristles will fade halfway approximately within 3 months, signaling it is time to change the brush head.

Replacement brush heads are available at your retailer or Oral-B Braun Service Centres.

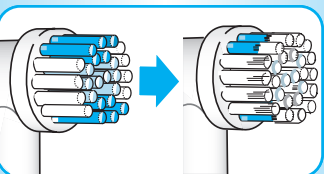After use

- Rinse your brush head thoroughly under running water for several seconds with the handle switched on.
- Switch off the handle and remove the brush head from the handle. Clean both parts separately under running water, then wipe them dry before you reassemble the brush and put the handle on the charger base.

Using your irrigator

Before using for the first time, fill the water container ⑭ and allow the water to run through completely by holding the irrigator handle ⑬ over a washbasin.

1. Attach one of the irrigator jets ⑧ to the handle until it locks into place. Fill the water container with luke warm water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch ⑩. When doing so, the handle switch ⑫ should be in the lower position = water stop.
3. Select the water pressure with the pressure control dial ⑪ (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch ⑫ up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.

The irrigator jet offers two different functions depending on the position of the jet tip:







Single jet:

Jet tip in position (a) to clean the teeth and the interdental spaces and remove resistant food particles.

Multi jet:

Jet tip in position (b) to massage the gums and keep them healthy.

First, turn off the water stream by sliding the handle switch down (to avoid unintended water spraying). The jet can easily be switched from one function to the other by pressing the jet

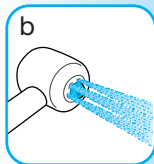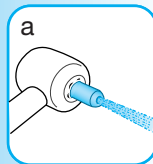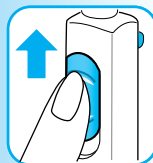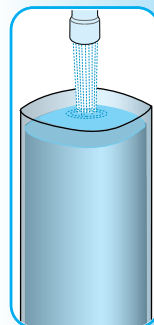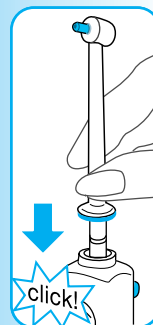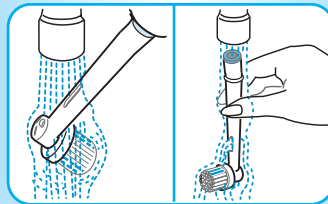
tip in or out. Begin with the single jet. Guide the water stream to the teeth and especially to the interdental spaces. Then massage the gums with the multi jet.

Max. operating time: 15 minutes

Cooling off time: 2 hours

To ensure the jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

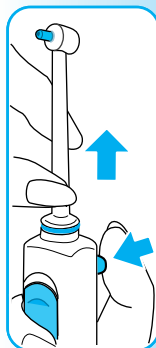After use

Always empty the water container entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended the use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator. To take off the irrigator jet, press the jet release button . The jets can be stored in the jet compartment. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth. The water container and the jet compartment  are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives EMC 2004/108/EC and Low Voltage 2006/95/EC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

For UK only:

This guarantee in no way affects your statutory rights.