
Oral-B®

ProfessionalCare

OxyJet-Center

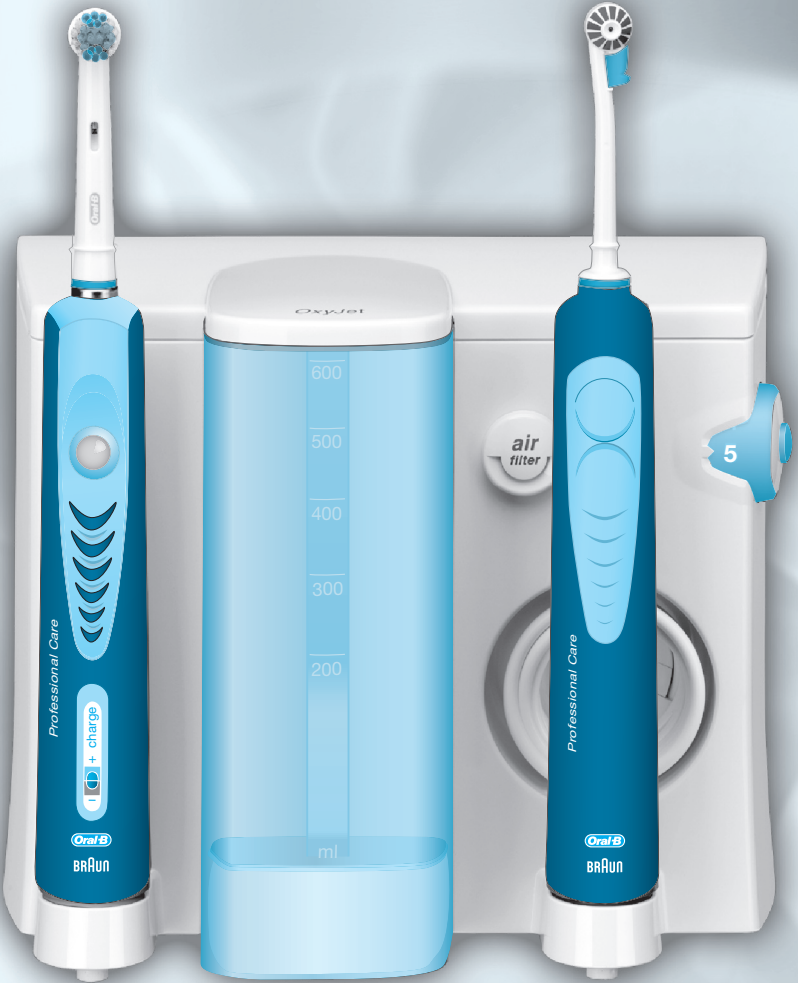powered
by **BRAUN**

Deutsch	4	DE AT	00 800 27 28 64 63
English	12	CH	08 44 - 88 40 10
Français	19	GB	0 800 731 1792
Español	27	IE	1 800 509 448
Português	35	FR	0 800 944 803 (appel gratuit depuis un poste fixe)
Italiano	43	BE	0 800 14 592
Nederlands	51	LU	4 00 50 51
Ελληνικά	59	ES	901 11 61 84
		PT	808 20 00 33
		IT	800 124 600
		NL	0800-4 45 53 88
		GR	800 801 3457

Internet:

www.oralb.com
www.braun.com
www.service.braun.com

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The **Oral-B ProfessionalCare OxyJet Center** is a dentist-inspired oral-care system for outstanding cleaning, gum care, whiter teeth and freshening.

The **Oral-B ProfessionalCare** toothbrush delivers outstanding performance:

- 3D oscillating-pulsating technology – for outstanding cleaning and gum care.
- Technologically advanced accessories and brush heads to clean and polish teeth as well as freshen your tongue and mouth.
- Brushing modes let you customize your brushing routine to address your specific needs.
- Premium handle design and ergonomics for greater comfort and control.

The **Oral-B OxyJet** irrigator mixes air into water in the form of micro-bubbles. It helps refresh your mouth. At the same time, the irrigator removes remaining food particles from between your teeth and massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

- The appliance must be mounted in a way that it cannot fall into water or drop to the floor. The plug should always be accessible.
- Periodically check the cord for damage. A damaged or non-functioning unit should no longer be used. If the cord is damaged, take the charging unit to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age 12.
- This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by a person responsible for their safety.
- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not place the appliance in water or other liquid.
- Do not reach for an appliance that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Braun Service Centre. Do not disassemble the product except when disposing of the battery.
- When taking out the battery from the toothbrush handle for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Do not insert any object into any opening of the appliance.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.

- Use this product only for its intended use as described in this manual. Do not use attachments or charges not recommended by the manufacturer.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

Toothbrush

- ① Charger base
- ② Charge indicator
- ③ Toothbrush handle (rechargeable)
- ④ On/off button
- ⑤ Brush head (see «Brush heads and accessories»)
- ⑥ Brush head compartment

Irrigator

- ⑦ Jet compartment
- ⑧ Irrigator jet
- ⑨ Jet switch
- ⑩ Jet release button
- ⑪ On/off switch
- ⑫ Water pressure control dial
- ⑬ Handle switch (water stop and go)
- ⑭ Irrigator handle
- ⑮ Water container (600 ml)
- ⑯ Cord storage

Specifications:


For voltage specifications please refer to the bottom of the charging unit.

Using your toothbrush

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and is designed for use in the bathroom.

- Plug the appliance into an electrical outlet and place the toothbrush handle ③ on the charger base ①.
- The **green** charge indicator ② light flashes while the handle is being charged.
- Once the handle is fully charged, the green charge indicator light shines continuously for 5 seconds and then it turns off. A full charge takes at least 18 hours and enables approximately one week of regular brushing (twice a day, 2 minutes). To check the charging status when the charging light is off, shortly lift up and return your handle on the charging unit. The charge indicator light will shine permanently for 5 seconds and turn off when your handle is fully charged.
- The **red** charge indicator ② light starts blinking when the rechargeable battery is low. It is now time to recharge your toothbrush.

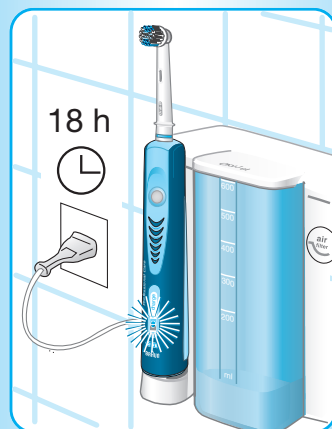OC 19.555.1

«Clean»

OC 19.575.3 X

«Clean»

«Sensitive»

«Polish»

- For everyday use, the toothbrush handle can be stored on the charger base to maintain it at full power. Overcharging is impossible. However, for environmental reasons, Braun Oral-B recommends to unplug the charging unit until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, refrain from placing the handle on the charger base and fully discharge it through regular use at least every 6 months.

Brushing modes

Your toothbrush offers different brushing modes to best meet your individual needs:


- «Clean» – Exceptional mouth cleaning.
- «Sensitive» – Gentle, yet thorough cleaning for sensitive areas.
- «Polish» – Exceptional polishing and natural whitening by removal of surface stains for occasional or everyday use.

How to switch between modes:

- Push the on/off button ④ once, and your toothbrush automatically starts in the «Clean» mode.
- To shift from the «Clean» mode to the «Sensitive» mode, push the on/off button once more.
- If your toothbrush has a «Polish» mode, pushing the on/off button again will bring you from «Sensitive» to «Polish» mode.
- To turn your brush off from the «Polish» mode, push the on/off button once more.
- If you wish to turn your toothbrush off in the «Clean» or «Sensitive» mode, push and hold down the on/off button until the motor stops.

Brushing technique

- Attach one of the brush heads ⑤ to the handle. They can be used with any kind of toothpaste.
- To avoid splashing, guide the brush head to your teeth before switching the appliance on.
- Guide the brush head slowly from tooth to tooth.
- Hold the brush head in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

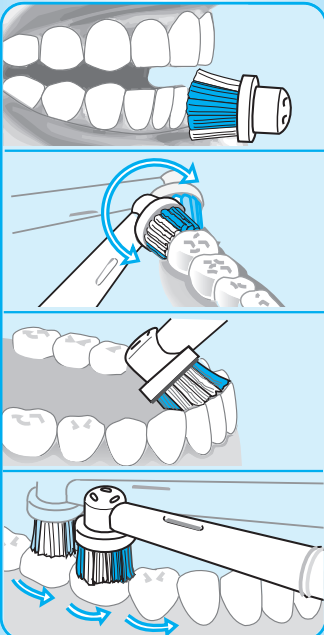Pressure Sensor

Oral-B ProfessionalCare has a pressure control feature to help you brush gently. This helps protect your gums from excessive brushing and helps make your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition to feeling it, you can also hear a different sound.

Professional Timer

A short stuttering sound at 30-second intervals, which you can hear and feel, reminds you to brush the four quadrants of your mouth equally for a thorough consistent clean.

2-Minute Timer

A long stuttering sound indicates when the dentist recommended 2-minute brushing time has been reached.

Both timers memorise the elapsed brushing time, even when the handle is briefly switched off during brushing. Only when pausing longer than 30 sec. the timer resets to the initial starting point.

In the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and gums, Oral-B recommends using the «Sensitive» mode (depending on model).

Brush heads and accessories




The **Oral-B ProfessionalCare** range offers you a variety of brush heads and accessories, allowing you to select those which best meet your individual needs. They are all designed for daily use:


I Oral-B PrecisionClean brush head

The Oral-B PrecisionClean brush head features outstanding oscillating-rotating action to effectively clean your teeth, clean along the gumline and in hard-to-reach areas such as behind teeth. The compact round brush head surrounds each tooth for a thorough clean. The green PrecisionClean bristles bend, enabling both gentle brushing and deep reach of the blue Interdental Tips between teeth. Apply the brushing technique described in the «Brushing technique» section.

II Oral-B Sensitive brush head (depending on model)


The Oral-B Sensitive brush head provides thorough cleaning for sensitive teeth and gums. It has an extra soft bristle texture that is gentle on teeth and gums. Apply the same brushing technique described in the «Brushing technique» section.

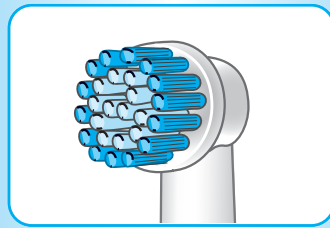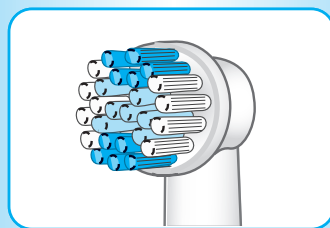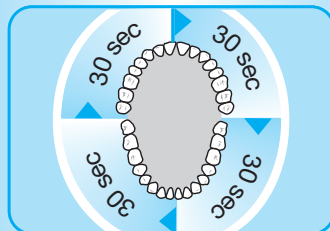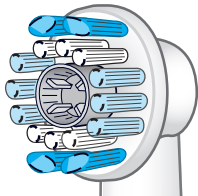
III Oral-B ProBright brush head (depending on model)

The Oral-B ProBright brush head is specifically designed to polish and naturally whiten your teeth by removing surface stains that can come from coffee, tea and tobacco. It features the oscillating-rotating head action of the Oral-B PrecisionClean brush head to effectively clean teeth, combined with the special polishing cup for natural brightening. Apply the same brushing technique described in the «Brushing technique» section. We do not recommend using the Oral-B ProBright brush head with braces. Children under 12 should not use the Oral-B ProBright brush head.



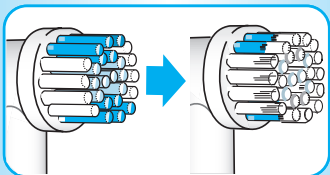
All three brush heads feature blue INDICATOR bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue colour will fade halfway approximately within 3 months, signaling it is time to change the brush head.

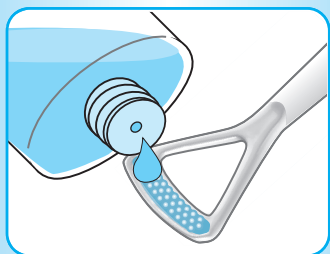IV Tongue Cleaner (depending on model)

The Tongue Cleaner helps freshen breath. It is designed to help remove oral malodour and can be used every day – before or after brushing. You can use mouthwash with the Tongue Cleaner. Attach it to the handle, turn the appliance on and drag it across the tongue from back to front. We recommend using the «Sensitive» mode for enhanced control during use. Rinse after use. Replace the Tongue Cleaner when the material becomes worn.

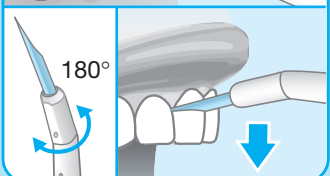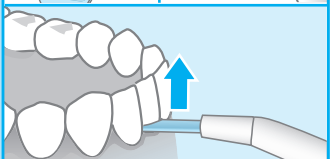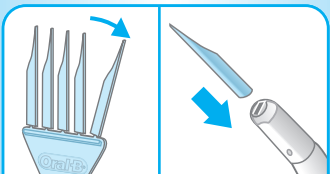V Interdental Cleaner (depending on model)

The Interdental Cleaner effectively cleans between the teeth and helps keep your gums healthy. It should be used before brushing. Break one pick cleaner off the base and insert it into the slot in the neck of the white pick holder. Attach the Interdental Cleaner to the handle and turn the appliance on. Depending on which part of your mouth you are cleaning, the pick holder should be adjusted by turning it, so that the blunt side always faces your gums. Guide the cleaner between your teeth and clean systematically from space to space. We recommend using the «Sensitive» mode for enhanced control during use. A pick cleaner can be used approximately two to three times until it becomes too soft to effectively fit and clean between the teeth. Replacement cleaners are available separately.

Oral-B ProfessionalCare replacement brush heads and accessories are available in stores that sell Oral-B products.

After use

- Rinse your brush head or accessory thoroughly under running water for several seconds with the handle switched on.
- Switch off the handle and remove the brush head or accessory from the handle. Clean both parts separately under running water, then wipe them dry before you reassemble the brush and put the handle on the charger base.

Using your irrigator

Before using for the first time, fill the water container ¹⁶ and allow the water to run through completely by holding the irrigator handle ¹⁴ over a washbasin.







1. Attach one of the irrigator jets ⁸ to the handle until it locks into place. Fill the water container with luke warm water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch ¹¹. When doing so, the handle switch ¹³ should be in the lower position = water stop.
3. Select the water pressure with the pressure control dial ¹² (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch ¹³ up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.



You can choose between two different jet functions by sliding the jet switch ⁹ down or up (before doing so, turn off the water stream by sliding the handle switch down):

Rotating stream for general cleaning and gum massage (jet switch in lower position): The rotation of the water stream and its micro-bubbles is designed for easy penetration.

Mono stream: straight water stream for site-specific cleaning (jet switch in upper position).

Max. operating time: 15 minutes
Cooling off time: 2 hours

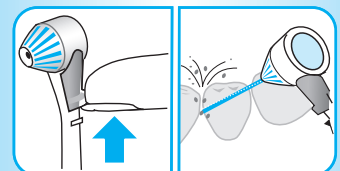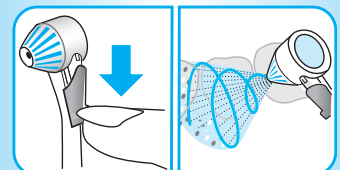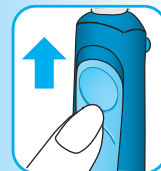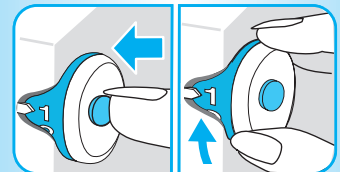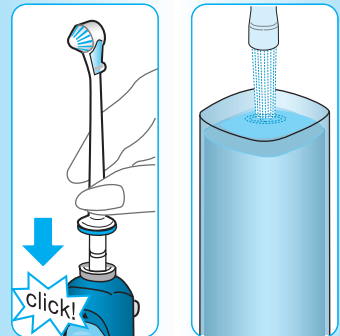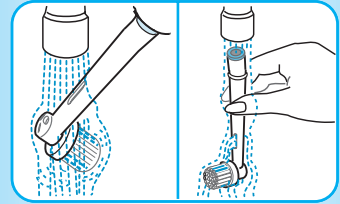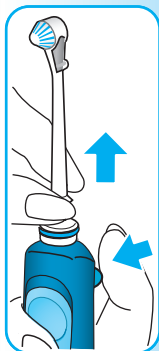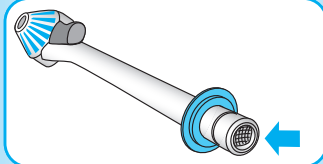To ensure the jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

After use

Always empty the water container entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended the use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch ⑩ to turn off the irrigator.

To take off the irrigator jet, press the jet release button ⑪. The jets can be stored in the jet compartment ⑦. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth. The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives EMC 2004/108/EC and Low Voltage 2006/95/EC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

For UK only:

This guarantee in no way affects your statutory rights.