Oral-B®

Professional Care

powered by **BRAUN**

Deutsch	4	(D) (A)	00 800 27 28 64 63
English	13	(CH)	08 44 - 88 40 10
Français	22	(UK)	0 800 416 434
Español	33	(IRL)	1 800 509 448
Português	43	(F)	0 810 365 855
Italiano	52	(B)	0 800 14 592
Nederlands	60	(L)	4 00 50 51
Ελληνικά	70	(ES)	901 11 61 84
		(PT)	808 20 00 33
		(I)	800 124 600
		(NL)	0800-4 45 53 88
		(GR)	800 801 3456
		(HK)	852-25249377 (Audio Supplies Company Ltd.)

Internet:

www.oralb.com

www.braun.com

www.service.braun.com

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Oral-B ProfessionalCare has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. However, when using electrical products, basic safety precautions should always be followed. Please read carefully before use.


Important

- Periodically check the cord for damage. A damaged or non-functioning appliance shall no longer be used. Take it to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are given supervision by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.

Description

- a Oral-B PrecisionClean brush head
- b On/off button
- c Handle
- d Charge indicator
- e Charging unit
- f Brush head compartment
- g Wall holder
- h Range of brush heads and accessories

Connecting and charging

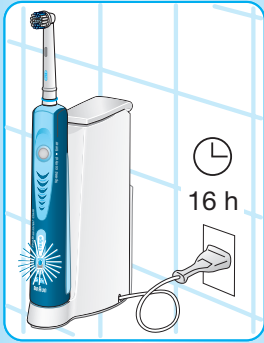
Your toothbrush has a waterproof handle, is electrically safe and is designed for use in the bathroom.

- Plug the charging unit (e) into an electrical outlet and place the toothbrush handle (c) on the charging unit.
- The green charge indicator (d) light flashes while the handle is being charged. If you then use your brush and place it back on the charging unit, the green charging light will again flash to let you know it is charging.
- When the handle is fully charged, the green charge indicator light shines continuously. A full charge takes at least 16 hours and enables approximately two weeks of regular brushing (twice a day, 2 minutes).
- The red charge indicator (d) light starts blinking when the rechargeable battery is low. Now it is time to recharge your toothbrush.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle by regular use at least every 6 months.

Specifications

Refer to the bottom of the charging unit for voltage specifications.

Using your toothbrush

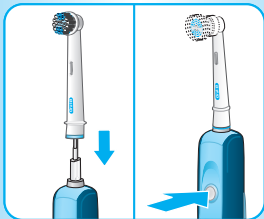Attach one of the brush heads (1–3) to the handle. They can be used with any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on.

Brushing modes

Your toothbrush offers different brushing modes to best meet your individual needs:

- «Clean» – Exceptional mouth cleaning.
- «Sensitive» – Gentle, yet thorough cleaning for sensitive areas.
- «Polish» – Exceptional polishing and natural whitening for occasional or everyday use.

How to switch between modes:

- Push the on/off button once, and your toothbrush automatically starts in the «Clean» mode.
- To shift from the «Clean» mode to the «Sensitive» mode, push the on/off button once more.
- If your toothbrush has a «Polish» mode, pushing the on/off button again will bring you from «Sensitive» to «Polish» mode.
- To turn your brush off from the «Polish» mode, push the on/off button once more.
- If you wish to turn your toothbrush off in the «Clean» or «Sensitive» mode, push and hold down the on/off button until the motor stops.

Outstanding Plaque Removal

For clean teeth and strong and healthy gums, **Oral-B ProfessionalCare** gently cleans in areas where plaque bacteria can collect. It provides an outstanding plaque removal even in hard-to-reach areas and is gentle on teeth and gums.

Brushing technique:

- Guide the brush head slowly from tooth to tooth.
- Hold the brush head in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

D 19.523.1

«Clean»

D 19.535.2

«Clean»


«Sensitive»

D 19.575.3 X

«Clean»


«Sensitive»

«Polish»

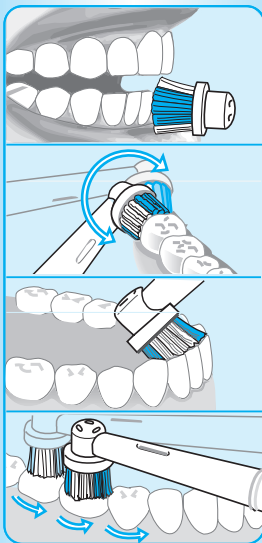Pressure Sensor

For optimal brushing results **Oral-B ProfessionalCare** has a pressure control feature to help you brush gently. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition to feeling it, you can also hear a different sound.

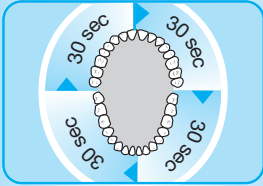Professional Timer

A short stuttering sound at 30-second intervals, which you can hear and feel, reminds you to brush the four quadrants of your mouth equally for a thorough consistent clean.

2-Minute Timer

A long stuttering sound indicates when the dentist recommended 2-minute brushing time has been reached.

Both timers memorise the elapsed brushing time, even when the handle is briefly switched off during brushing. Only when pausing longer than 30 sec. the timer resets to the initial starting point.


In the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or hygienist. If you have sensitive teeth and gums, Oral-B recommends using the «Sensitive» mode.

Individual brush heads and accessories

The **Oral-B ProfessionalCare** range offers you a variety of brush heads and accessories, allowing you to select those which best meet your individual needs. They are all designed for daily use:

1 Oral-B PrecisionClean brush head


The Oral-B PrecisionClean brush head features outstanding oscillating-rotating action to effectively clean your teeth, clean along the gumline and in hard-to-reach areas such as behind teeth. The compact round brush head surrounds each tooth for a thorough clean. The green PrecisionClean bristles bend, enabling both gentle brushing and deep reach of the blue Interdental Tips between teeth. Apply the brushing technique described in the «Outstanding Plaque Removal» section.

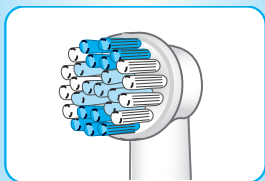2 Oral-B Sensitive brush head

(not with all models)


The Oral-B Sensitive brush head provides thorough cleaning for sensitive teeth and gums. It has an extra soft bristle texture that is gentle on teeth and gums. Apply the same brushing technique described in the «Outstanding Plaque Removal» section.

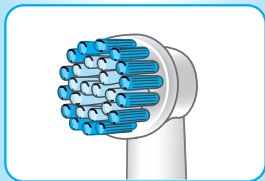
3 Oral-B ProBright brush head


(not with all models)

The Oral-B ProBright brush head is specifically designed to polish and naturally whiten your teeth by removing stains that can come from coffee, tea and tobacco. It features the oscillating-rotating head action of the Oral-B PrecisionClean brush head to effectively clean teeth, combined with the special polishing cup for natural brightening. Apply the same brushing technique described in the «Outstanding Plaque Removal» section. We do not recommend using the Oral-B ProBright brush head with braces. Children under 12 should not use the Oral-B ProBright brush head.

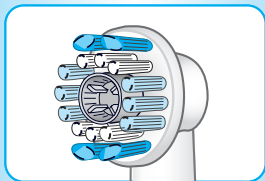All three brush heads feature blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue colour will fade halfway approximately within 3 months, signaling it is time to change the brush head.

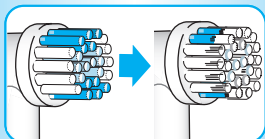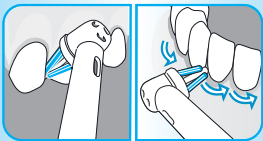
4 Interspace brush head

(not with all models)

The Interspace brush head is designed for cleaning between teeth as well as around bridges, crowns and implants.


To use, position brush head at a 45 degree angle to the gumline. Slowly guide the brush head along the gumline, pausing between teeth for excellent results.

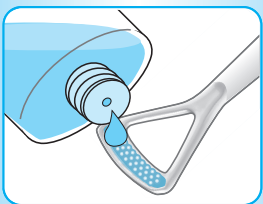5 Tongue Cleaner

(not with all models)

The Tongue Cleaner helps freshen breath. It is designed to help remove oral malodour and can be used every day – before or after brushing. You can use mouthwash with the Tongue Cleaner. Attach it to the handle, turn the appliance on and drag it across the tongue from back to front. We recommend using the «Sensitive» mode for enhanced control during use. Rinse after use. Replace the Tongue Cleaner when the material becomes worn.

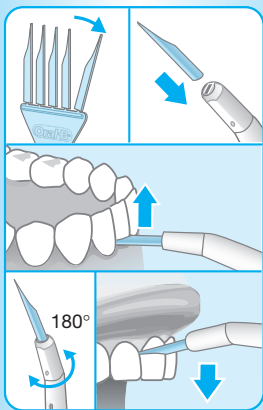6 Interdental Cleaner


(not with all models)

The Interdental Cleaner effectively cleans between the teeth and helps keep your gums healthy. It should be used before brushing. Break one pick cleaner off the base (7) and insert it into the slot in the neck of the white pick holder. Attach the Interdental Cleaner to the handle and turn the appliance on. Depending on which part of your mouth you are cleaning, the pick holder should be adjusted by turning it, so that the blunt side always faces your gums. Guide the cleaner between your teeth and clean systematically from space to space. We recommend using the «Sensitive» mode for enhanced control during use. A pick cleaner can be used approximately two to three times until it becomes too soft to effectively fit and clean between the teeth. Replacement cleaners are available separately.

Oral-B ProfessionalCare replacement brush heads and accessories are available in stores that sell Oral-B products.

Cleaning recommendations

- After brushing, rinse your brush head or accessory thoroughly under running water for several seconds with the handle switched on.
- Switch off the handle and remove the brush head or accessory from the handle. Clean both parts separately under running water, then wipe them dry before you reassemble the brush and put the handle on the charging unit.
- From time to time, also clean the charging unit and the brush head compartment (D 19.575.3 X, D 19.535.2 models only) with a damp cloth.
- For cleaning purposes, the brush head compartment (D 19.575.3 X, D 19.535.2 models only) is removable.

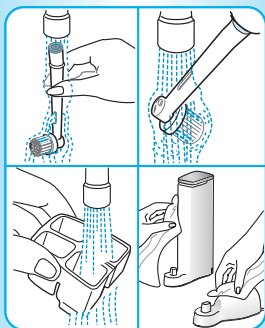Subject to change without notice.

This product conforms to the European Directives EMC 2004/108/EC and Low Voltage 2006/95/EC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

For UK only:

This guarantee in no way affects your statutory rights.

The 30-Day Money Back Guarantee

Take the challenge! Experience the **Oral-B ProfessionalCare** difference! Satisfaction guaranteed or your money back.

Oral-B ProfessionalCare provides outstanding cleaning, natural whitening by removing stains and takes care of your gums. In fact, we are so convinced you'll love your **Oral-B ProfessionalCare** that we give you this special offer.

Try **Oral-B ProfessionalCare** for 30 days starting from the day of purchase. If you are not 100% satisfied, receive a full refund.

For a full refund, please send us your handle, refills, and charging unit along with your name, address, contact phone number and the original till receipt indicating purchase price and date.

United Kingdom / Ireland:
«Oral-B 30-Day Money Back
Guarantee»
Dept PX415,
PO Box 10,
Diss, Norfolk,
IP98 1SP

For further details, please contact our helpline (see page 2) to be referred to the authorised Oral-B Braun Service Centre closest to you.