CUPS EACH TOOTH
FOR A
3D CLEAN

Oral-B®

TRIUMPH® powered by **BRAUN**

4000

• PRESSURE SENSOR

Deutsch	4
English	11
Français	18
Español	26
Português	33
Italiano	40
Nederlands	47
Dansk / Norsk	54
Svenska	61
Suomi	68
Ελληνικά	75

(D) (A)	00 800 27 28 64 63
(CH)	08 44 - 88 40 10
(UK)	0 800 731 1792
(IRL)	1 800 509 448
(F)	0 810 365 855
(BE)	0 800 14 592
(L)	4 00 50 51
(ES)	901 11 61 84
(PT)	808 20 00 33
(I)	800 124 600
(NL)	0800-4 45 53 88
(DK)	70 15 00 13
(N)	22-63-00-93
(S)	020 - 21 33 21
(FIN)	020 377 877
(GR)	800 801 3456
(HK)	852-25249377 (Audio Supplies Company Ltd.)

Internet:

www.oralb.com
www.braun.com
www.service.braun.com

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Oral-B Triumph® has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. When using electrical products, however, some basic safety precautions should always be followed.

IMPORTANT

- The Oral-B Triumph is provided with a SmartPlug (integrated Safety Extra Low Voltage) power supply. Do not exchange or manipulate any part of it, as there is a risk of electric shock. Periodically check the SmartPlug for damage and if the appliance is damaged or non-functional, stop use of the product. Take damaged units to an Oral-B Braun® Service Centre.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are supervised by a person responsible for their safety.
- If the product is dropped, the brush head should be replaced before the next use, even if no damage is visible.
- Keep the SmartPlug away from water or other liquid.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Braun Service Centre. Do not disassemble the product, except when disposing of the battery.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Do not touch the SmartPlug with wet hands. This can cause electric shock.
- Use this product only for its intended use as described in this manual. Do not use attachments or chargers not recommended.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description



- a Oral-B FlossAction® brush head
- b On/off button
- c Brushing mode button
- d Handle
- e Handle display
- f SmartPlug charger (includes charging ring and SmartPlug)
- g Base plate
- h Protective cover
- i Oral-B® ProBright brush head (not with all models)

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

To assemble the charging unit, insert the protective cover (h) on the back of the base plate (g). Then, place the base plate on the charging ring (f) until it clicks into place.

- Plug the SmartPlug into an electrical outlet and place the toothbrush handle (d) on the charging unit.

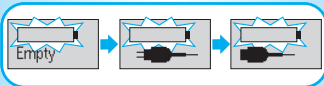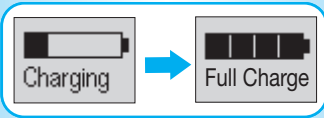- A battery symbol should appear in the handle display (e). The battery symbol will fill up while the handle is being charged and «Charging» will be displayed. When the handle is fully charged, the handle display shows «Full Charge». A full charge takes at least 10 hours and enables approximately 10 days of regular brushing (twice a day, 2 minutes).

Please note: If the battery symbol is not displayed immediately, continue charging. It will appear in 10–15 minutes.

- When the rechargeable battery is low, the bars of the battery symbol will disappear and «Low» will be displayed. Your toothbrush will slow down and then stop. It is time to recharge Oral-B Triumph. When the battery is fully discharged, «Empty» is displayed. It will then need a 15-minute charge before you can use it for a 2-minute brushing.

- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. There is no risk of overcharging the battery.

- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

For electric specifications, see printing on the SmartPlug.

Personalizing your toothbrush

You can program your toothbrush so that it meets your specific needs. You can choose the language and the timer that works best for you.

Select your language

Your toothbrush comes with the English language activated. If you wish to switch to another language, proceed as follows: Press and hold the brushing mode button (c) until «English» is displayed. You can scroll through 13 different languages by briefly pressing the on/off button (b). To confirm the desired language, briefly press the brushing mode button again. This will also lead you to the timer choice menu.

Select your timer

You can choose between the «Professional» or «2-Minute» timer to meet your specific needs.

The «2-Minute Timer» indicates by a long stuttering sound when the dentist recommended 2-minute brushing time has been reached.

The «Professional Timer» helps you brush all four quadrants of your mouth equally. It signals with a short stuttering sound at 30-second intervals (45 seconds intervals within «Deep Clean» mode) that it is time to move to the next quadrant of your mouth, followed by a long stuttering sound when you reach the recommended two minute mark.



The total brushing time can be extended to 3 minutes by using the «Deep Clean» mode (see «Brushing modes»).


Both timers memorize the elapsed brushing time, even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds the timer resets to the initial starting point.

Your toothbrush comes with the «Professional Timer» activated. If you wish to switch to the «2-Minute Timer» press the on/off button (b) after you have confirmed the language. Then, press the brushing mode button (c) to confirm your choice.

Brushing

You can start your toothbrush by pressing the on/off button (b). To avoid splashing, guide the brush head to your teeth before switching the appliance on.

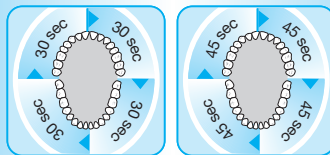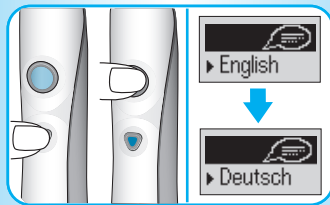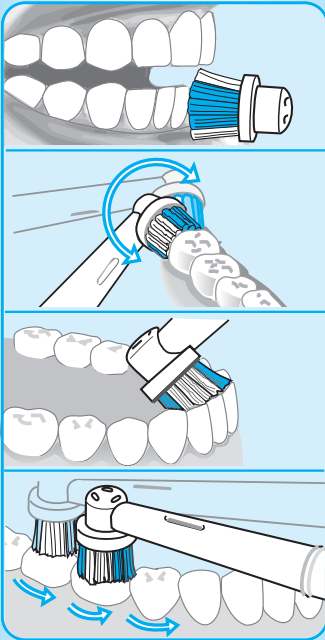Brushing technique

- Guide the brush head slowly from tooth to tooth.
- Hold the brush head in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

Oral-B Triumph has a pressure control feature to help ensure that the right degree of pressure is applied. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition to feeling the change, you can also hear a different sound.

In the first few days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or hygienist.


Brushing modes (not with all models)

Your toothbrush offers different brushing modes for varying oral care needs:

- «Daily Clean» – Exceptional mouth cleaning.
- «Sensitive» – Gentle, yet thorough cleaning for sensitive areas.
- «Massage» – Gentle massage of gums.
- «Polish» – Polishing for occasional or everyday use.
- «Deep Clean» – Exceptional mouth cleaning with extended brushing time. Within this mode the brushing time has been extended to 45 seconds per quadrant for a total brushing time of 3 minutes.

How to switch between brushing modes:

- When pressing the on/off button your toothbrush automatically starts in the «Daily Clean» mode.
- To switch to the other modes, successively press the mode button (c).
- To return from any mode to the «Daily Clean» mode, press and hold the mode button.
- To turn your brush off, press the on/off button.

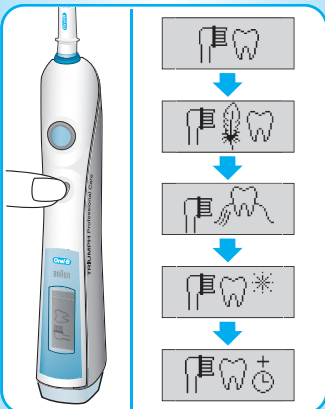During brushing, the handle display shows the brushing mode you have selected, followed by your actual brushing time. When finished brushing, the total brushing time will be displayed for 30 seconds before the display switches back to the battery charge level. To directly get the battery charge level display, you can clear the brushing time display by pushing the brushing mode button (c).

Please note: The «Polish» mode is best used with an Oral-B ProBright (i) brush head (not with all models).

Brush heads

The Oral-B FlossAction brush head (a) is recommended for everyday deep cleaning. It features micropulse bristles that get deep in hard-to-reach areas between your teeth, gently removing hidden plaque.





For occasional or everyday polishing, please use the Oral-B ProBright (i) brush head (not with all models). It features a dentist-inspired polishing cup to remove surface stains from coffee, tea and tobacco for effective brightening. Please note that children under 12 years old should not use the Oral-B ProBright brush head.


Do not use the Oral-B FlossAction or Oral-B ProBright brush head with braces. We recommend using the Oral-B Ortho brush head, which is specifically designed to clean around braces and wires.

All brush heads feature blue Indicator® bristles helping you to monitor brush head replacement need. With proper brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, signalling that it is time to change the brush head.


Cleaning recommendations

- After brushing, rinse your brush head thoroughly under running water.
- Remove the brush head from the handle and clean both parts separately under running water. Then, wipe them dry before you put the handle on the charging unit.

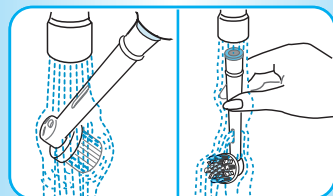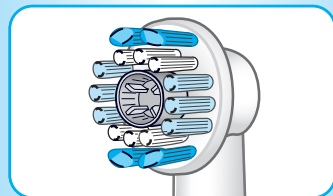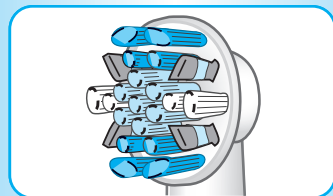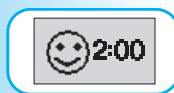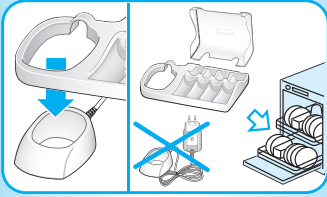- Before cleaning your charging unit, make sure you disassemble it first. **The charging ring and the SmartPlug should never be placed in the dishwasher.** However, all other parts (base plate and protective cover) are dishwasher-safe.

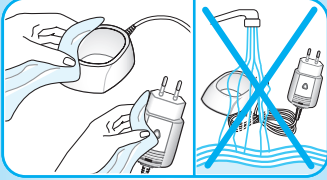- From time to time, you can also clean the charging ring and SmartPlug unit with a damp cloth. The SmartPlug charger should never be placed in water.

Subject to change without notice.

This product conforms to the European Directives EMC 2004/108/EC and Low Voltage 2006/95/EC.

Environmental notice

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at approved recycling or disposing locations in accordance with local regulations.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.

Satisfaction Guaranteed!

Take the challenge! Experience the Oral-B Triumph difference!

Oral-B Triumph is Oral-B's technologically advanced rechargeable toothbrush for most outstanding cleaning, polishing, gentle gum care and freshening of your mouth. We are so convinced you'll fall in love with your new toothbrush that we give you this risk-free offer.

Try Oral-B Triumph for 30 days, starting from the day of purchase. If you are not 100% satisfied, receive a full refund. Guaranteed. For a full refund, please return your handle, refills and charging unit with your original cash register receipt to the address below within 60 days of purchase.

Great Britain / Ireland:
«Oral-B 30-Day Challenge»
Braun Cashback (BC 024)
PO Box 504
Leicester
LE94 0AE

Terms and conditions apply, for further details, please contact our customer service 0800 7311 792 to be referred to the authorized Oral-B Braun Service Centre closest to you.

This does not affect your statutory rights.