Oral-B®

TRIUMPH® powered by BRAUN

4000

Deutsch	4	DE AT	00 800 27 28 64 63
English	7	CH	08 44 - 88 40 10
Français	10	GB	0 800 731 1792
Español	14	IE	1 800 509 448
Português	17	FR	0 800 944 803 (appel gratuit depuis un poste fixe)
Italiano	20		
Nederlands	23	BE	0 800 14 592
Dansk / Norsk	26	LU	4 00 50 51
Svenska	29	ES	901 11 61 84
Suomi	32	PT	808 20 00 33
Ελληνικά	35	IT	800 124 600
		NL	0800-4 45 53 88
		DK	70 15 00 13
		NO	22 63 00 93
		SE	020 - 21 33 21
		FI	020 377 877
		GR	800 801 3457

Internet:

www.oralb.com
www.braun.com
www.service.braun.com

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English

Welcome to Oral-B!

Before operating this toothbrush, please read the instructions and save this manual for future reference.

IMPORTANT

- Periodically check the cord for damage. If the cord is damaged, take the charging unit to an Oral-B Braun Service Centre. A damaged or non-functioning unit must no longer be used.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.
- Children should be supervised to ensure that they do not play with the appliance.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service Centre. Do not disassemble the product except when disposing of the battery.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Do not insert any object into any opening of the appliance/charger.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

- a Brush head
- b Pressure sensor light
- c On/off button
- d Brushing mode button
- e Handle

- f Brushing mode display
- g Charge indicator
- h Charging unit (basic charger and brush head compartment with protective cover)

Connecting and charging

- Plug the charger into an electrical outlet. For electric specifications, see printing on the charger.
- The flashing green lights indicate that the toothbrush is charging; once it is fully charged, the lights turn off. A full charge enables approx. 10 days of regular brushing (twice a day, 2 minutes). Note: If the green lights do not flash immediately, continue charging. They will appear in 10–15 minutes.
- If the rechargeable battery is low the red low charge indicator light is flashing for some seconds when turning your toothbrush on/off. If the rechargeable battery is empty the motor will stop; it will need a 40-minute charge before you can use it for a 2-minute brushing.
- You can check the charging status for some seconds by turning your toothbrush on/off or placing it on the charging unit.

Using your toothbrush

Brushing technique

Wet brush head and apply toothpaste; you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on (c). Guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface. Start brushing the outsides, and then the insides and finally the chewing surfaces (picture 2, page 3). Brush all four quadrants of your mouth equally.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and gums, Oral-B recommends using the «Sensitive» mode optional in combination with Oral-B «Sensitive» brush head.

Professional Timer

The «Professional Timer» signals with short beep sounds at 30-seconds intervals to move to the next quadrant of your mouth (45 seconds intervals within «Deep Clean» mode) (picture 2, page 3). A single beep indicates completion of the first quadrant, two beeps the second one, three beeps the third one and four beeps sound when you reach the recommended 2-minute brushing time.

The «Professional Timer» memorizes the elapsed brushing time, even when the handle is briefly switched off during brushing. The timer resets if pausing longer than 30 seconds or briefly pressing the mode button (d) during pause.

Pressure sensor

If you put excessive pressure on your teeth, the red pressure sensor light (b) will light up (picture 2, page 3). Additionally, when brushing in «Daily Clean» or «Deep Clean» mode, the brushing mode will change to the «Sensitive» mode. Periodically check the operation of the pressure sensor by pressing moderately on the brush head during use.

Brushing modes

Your toothbrush offers different brushing modes:


- «Daily Clean» – Standard mode for daily mouth cleaning
- «Deep Clean» – Mouth cleaning with extended brushing time of 45 seconds per quadrant (3 minute total brushing time)
- «Sensitive» – Gentle, yet thorough cleaning for sensitive areas
- «Polish/Whitening» – Polishing for occasional or everyday use

Your toothbrush automatically starts in the «Daily Clean» mode. To switch to the other modes, successively press the mode button (picture 3, page 3). To return from any mode back to the «Daily Clean» mode, press and hold the mode button. The brushing mode is memorized even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds or briefly pressing the mode button (d) during pause the brushing mode resets to the «Daily Clean» mode.

Brush heads


Oral-B offers you a variety of different brush heads that fit your Oral-B Triumph toothbrush.

Oral-B FlossAction brush head

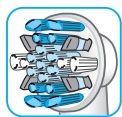features micropulse bristles that allow for an outstanding plaque removal in the interdental areas.

Oral-B Precision Clean

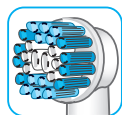
cups each tooth with a curved bristle trim and reaches deep between teeth at the same time.

Oral-B Sensitive brush head

features an extra soft bristle structure that is gentle on teeth and gums.

Oral-B 3DWhite brush head

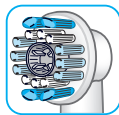features a special polishing cup to naturally whiten your teeth. Please note that children under 12 years old should not use the Oral-B 3DWhite brush head.

All brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway down approximately in a 3-month period. If the bristles splay before the colour recedes, you may be exerting too much pressure on teeth and gums.

We do not recommend using the Oral-B FlossAction or Oral-B 3DWhite brush head with braces. You may use the Oral-B Ortho brush head, specifically designed to clean around braces and wires.

Cleaning recommendations

After brushing, rinse your brush head under running water with the handle switched on. Switch off and remove the brush head. Clean handle and brush head separately; wipe them dry before reassembling. Disassemble the charging unit before cleaning. Brush head compartment and protective cover are dishwasher-safe. Clean the basic charger with a damp cloth only (picture 4, page 3).

Subject to change without notice.

Environmental Notice

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at approved recycling or disposing locations in accordance with local regulations.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.

Satisfaction Guaranteed or your money back!

Take the 30-day Challenge!

Oral-B Triumph is Oral-B's most technologically advanced rechargeable toothbrush for outstanding cleaning and whitening, gentle gum care and freshening of your mouth. We are so convinced you'll fall in love with your new toothbrush that we give you this risk-free offer.

Try Oral-B Triumph for 30 days, starting from the day of purchase. If you are not 100% satisfied, receive a full refund. No questions asked. For a full refund, please return your handle, refills and charging unit with your original cash register receipt to the address below within 30 days of purchase.

Great Britain / Ireland:

«Oral-B 30-Day Challenge»
Braun Cashback (BC 024)
PO Box 504
Leicester
LE94 0AE

UK warranty:

FREEPOST RSHJ-KRSA-UHGE
Braun Product Service
Olympic Shaver Centre Ltd
PO BOX 9799
Calverton
Nottingham
NG14 6WN

ROI warranty:

Braun Product Service
Olympic Shaver Centre Ltd
Dept AA1758
PO BOX 6151
FREEPOST
Dublin 2

For further details, please contact our customer service 0800 7311 792 to be referred to the authorized Oral-B Braun Service Centre closest to you.