Oral-B®

TRIUMPH® powered by BRAUN

5000

with Wireless Smart Guide

Deutsch	6	DE AT	00 800 27 28 64 63
English	12	CH	08 44 - 88 40 10
Français	17	GB	0 800 731 1792
Español	23	IE	1 800 509 448
Português	28	FR	0 800 944 803 (appel gratuit depuis un poste fixe)
Italiano	34		
Nederlands	39	BE	0 800 14 592
Dansk / Norsk	44	LU	800 21172
Svenska	50	ES	901 11 61 84
Suomi	55	PT	808 20 00 33
Ελληνικά	60	IT	800 124 600
		NL	0800-4 45 53 88
		DK	70 15 00 13
		NO	22 63 00 93
		SE	020 - 21 33 21
		FI	020 377 877
		GR	800 801 3457

Internet:

www.oralb.com
www.braun.com
www.service.braun.com

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English

Welcome to Oral-B!

Before operating this toothbrush, please read these instructions and save this manual for future reference.

IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- This product is not intended for use by children under age three. Children from 3 to 14 years, persons with reduced physical, sensory or mental capabilities, or persons with lack of experience or knowledge of the product can use the toothbrush, but should have supervision and/or instruction concerning use of the appliance in a safe way, and understand the hazards involved.
- Cleaning and maintenance should not be performed by children.
- Children should not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

SmartGuide batteries

- Batteries may leak if empty or not used for a long time. Avoid skin contact whilst handling leaking batteries.
- Keep away from children. Do not swallow. If swallowed, seek medical advice at once.
- Remove fully discharged batteries immediately. Replace all batteries at the same time. Do not mix used and new batteries, different brands or types. Do not disassemble, recharge or dispose of in fire.

WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water, liquid or store it where it can fall or be pulled into a tub or sink. Do not reach for it if it has fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Important Information

- Your Oral-B Triumph uses an interactive wireless display. To avoid electromagnetic interference and/or compatibility conflicts, deactivate the radio transmission before using it in restricted surroundings, such as aircrafts or specially marked areas in hospitals. Deactivate radio transmission by simultaneously pressing the on/off and mode button for 3 seconds until the radio transmission display (f) goes off. Follow the same procedure to reactivate transmission again.
- Persons with pacemakers should always keep the toothbrush more than six inches (15 centimeters) from the pacemaker while turned on. Anytime you suspect interference is taking place, deactivate the radio transmission of your toothbrush.

Description

- a Brush head
- b Pressure sensor light
- c On/off button
- d Brushing mode button
- e Handle
- f Radio transmission display
- g Charge indicator
- h Charging unit (basic charger and brush head compartment with protective cover)
- i SmartGuide
- j Wall mount
- k Adhesive tape

Connecting and charging

- Plug the charger into an electrical outlet. For electric specifications, see printing on the charger.

- The flashing green lights indicate that the toothbrush is charging (picture 1, page 5); once it is fully charged, the lights turn off. A full charge enables up to 10 days of regular brushing (twice a day, 2 minutes).
Note: If the green lights do not flash immediately, continue charging. They will appear in 10–15 minutes.
- If the rechargeable battery is low the red low charge indicator light is flashing for some seconds when turning your toothbrush on/off. If the rechargeable battery is empty the motor will stop; it will need a 40-minute charge before you can use it for a 2-minute brushing.
- You can check the charging status for some seconds by turning your toothbrush on/off or placing it on the charging unit.

SmartGuide adjustments

The SmartGuide (i) comes with a demonstration mode activated. Deactivate this mode by removing the battery compartment cover in the back (picture A, page 4) and push the «set» or «h/min» button. A flashing «12:00» indicates that you can now adjust the clock time.

Clock Time

Press and hold the «set» button until «12h» is displayed. By briefly pressing the «h/min» button, you can toggle between «12h» or «24h» daytime format. Push «set» to confirm your choice. Then the hour digit starts flashing. Press the «h/min» button until the right hour appears and confirm with «set». Proceed same way to set both minute digits (picture B, page 4). If you prefer not using the clock feature, the flashing «12:00» will disappear automatically after 24h.

SmartGuide Icons

Displayed Icon Meaning

«Daily Clean» mode

«Deep Clean» mode

«Sensitive» mode

«Polish» mode

«Massage» mode

Too much pressure is applied.
Reduce the brushing force.

«Quadrant circle»
Move to the next quadrant of your mouth when the relevant segment flashes.

You reached the dentist recommended 2-minute brushing time (3-minute brushing time within the Deep Clean mode).

You will receive a star with each completed quadrant.

Smart Guide Batteries are low.
Change batteries of the SmartGuide.

Using your toothbrush

Brushing technique

Wet brush head and apply toothpaste; you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on.

When brushing your teeth with one of our rotating brush heads (find out more in chapter Brush heads, page 14), guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface (picture 3, page 5).

When using our TriZone brush heads (find out more in chapter Brush heads, page 14) place the toothbrush bristles against the teeth at a slight angle towards the gumline. Apply light pressure and start

brushing in back and forth motions, just like you would do with a manual toothbrush. With all brush heads start brushing the outsides, then the insides and finally the chewing surfaces. Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and/or gums, Oral-B recommends using the «Sensitive» mode in combination with an Oral-B Sensitive brush head.

Timer

You can choose between «2-Minute» or «Professional» timer and «Count up» or «Count down» timer.

The «2-Minute» timer signals with four beeps and a Smiley on the SmartGuide that the dentist recommended 2-minute brushing time has been reached.

The «Professional» timer signals with short beep sounds at 30-seconds intervals to move to the next quadrant of your mouth (45 seconds intervals within «Deep Clean» mode) (picture 3, page 5). A single beep indicates completion of the first quadrant, two beeps the second one, three beeps the third one and four beeps sound when you reach the recommended 2-minute brushing time.

The timer memorizes the elapsed brushing time, even when the handle is briefly switched off during brushing. The timer resets if pausing longer than 30 seconds or briefly pressing the mode button (d) during pause.

The «Count Up» timer displays your actual brushing time. The «Count down» timer indicates the remaining brushing time.

Personalizing the timer

Your toothbrush comes with the «Professional» and «Count up» timer activated. To change this, press and hold the brushing mode button for 3 seconds until the SmartGuide shows «:30» and «timer». Briefly press the on/off button (c) to switch to the «2-Minute» timer. Push the mode button to confirm. The SmartGuide now shows «up» and «timer». Press the on/off button to switch to the «Count down» timer and confirm with the mode button. Leave the setup menu by pressing the mode button.

Pressure sensor

If you put excessive pressure on your teeth, the red pressure sensor light (b) will light up (picture 3, page 5). Additionally, when brushing in «Daily

Clean» or «Deep Clean» mode, the brushing mode will change to the «Sensitive» mode. Periodically check the operation of the pressure sensor by pressing moderately on the brush head during use.

Brushing modes

Your toothbrush offers different brushing modes:


- «Daily Clean» – Standard mode for daily mouth cleaning
- «Deep Clean» – Mouth cleaning with extended brushing time of 45 seconds per quadrant (3 minute total brushing time)
- «Sensitive» – Gentle, yet thorough cleaning for sensitive areas
- «Polish» – Polishing for occasional or everyday use
- «Massage» – Gentle massage of gums

Your toothbrush automatically starts in the «Daily Clean» mode. To switch to the other modes, successively press the mode button (picture 4, page 5). To return from any mode back to the «Daily Clean» mode, press and hold the mode button. The brushing mode is memorized even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds or briefly pressing the mode button (d) during pause the brushing mode resets to the «Daily Clean» mode.

Brush heads


Oral-B offers you a variety of different Oral-B brush heads that fit your Oral-B toothbrush handle. Our rotating brush heads can be used for tooth-by-tooth precision cleaning.

Oral-B FlossAction brush head

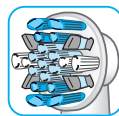features micropulse bristles that allow for an outstanding plaque removal in the interdental areas.

Oral-B Precision Clean brush head

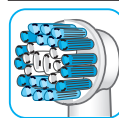
cup each tooth with a curved bristle trim and reaches deep between teeth at the same time.

Oral-B Sensitive brush head

features an extra soft bristle structure that is gentle on teeth and gums.

Oral-B 3DWhite brush head

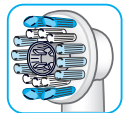features a special polishing cup to naturally whiten your teeth. Please note that children under 12 years old should not use the Oral-B 3DWhite brush head.

Oral-B PowerTip brush head

offers a deep clean in hard to reach areas, e.g. around bridges, crowns or implants.


Oral-B Ortho brush head

specifically designed for effective, convenient cleaning around orthodontic appliances.

Our dynamic TriZone brush head can be used with all common brushing techniques.

Oral-B TriZone brush head

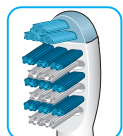Is designed to be used like a manual toothbrush while delivering outstanding plaque removal thanks to its triple zone cleaning action.

Most of our Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue colour will fade halfway down approximately in a 3-month period indicating that need to exchange your brush head to maintain a high level of plaque removal. If the bristles splay before the colour recedes, you may exert too much pressure on teeth and gums.

We do not recommend using the Oral-B FlossAction or Oral-B 3DWhite brush head with braces. You may use the Oral-B Ortho brush head, specifically designed to clean around braces and wires.

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.

- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

Synchronizing your toothbrush

To avoid interference of display messages when using more than one Oral-B Triumph toothbrush with SmartGuide in the same household, you need to assign the handles to their respective SmartGuide.

You can synchronize a maximum of 2 handles with one SmartGuide. To do so, push the «h/min» button in the battery compartment for at least 3 seconds. The SmartGuide shows «L---». Turn on your handle to start the synchronization process. It is finished when «L-1-» is displayed. Turn off the handle. To assign a second handle to the same SmartGuide turn on the other handle until «L-2-» is displayed. After having synchronized the second handle, «Full» will be displayed. Exit the synchronization mode by pressing the «h/min» or «set» button, otherwise it will end automatically after 30 seconds.

Cleaning recommendations

After brushing, rinse your brush head under running water with the handle switched on. Switch off and remove the brush head. Clean handle and brush head separately; wipe them dry before reassembling. Disassemble the charging unit before cleaning. Brush head compartment and protective cover are dishwasher-safe. Both basic charger and SmartGuide should be cleaned with a damp cloth only (picture 5, page 5).

Wall mount

Use the adhesive tape (k) if you prefer to mount the SmartGuide on a wall, mirror or other surface. Make sure the chosen surface is clean and dry. Position the wall mount (j) with the Oral-B logo readable and press it evenly and firmly to the wall. Wait 24h before sliding the SmartGuide into the fixed wall mount (picture C, page 4). The adhesive tape will not work on dirt-repellent surfaces.

Subject to change without notice.

This appliance complies with the requirements covered by the European Directive 1999/5/EC. Receiver Class 3, Duty Cycle Class 3. To find the complete EC declaration of conformity, please visit www.oralb.com/ce.

Environmental Notice

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Center or at approved recycling or disposing locations in accordance with local regulations.

Satisfaction Guaranteed or your money back!

Take the 30-day-challenge!

Oral-B Triumph with SmartGuide is Oral-B's most technologically advanced rechargeable toothbrush for outstanding cleaning and whitening, gentle gum care and freshening of your mouth. We are so convinced you'll fall in love with your new toothbrush that we give you this risk-free offer.

Try Oral-B Triumph with SmartGuide for 30 days, starting from the day of purchase. If you are not 100% satisfied, receive a full refund. No questions asked. For a full refund, please return your handle, refills and charging unit with your original cash register receipt to the address below within 30 days of purchase.

Oral B Money Back Guarantee
Dept. PX415
PO Box 10
Diss
Norfolk
IP98 3SP

For further details, please contact our customer service 0800 7311 792 to be referred to the authorised Oral-B Braun Service Centre closest to you.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

Please also read our additional warranty statements: «Warranty replacement brush heads» on page 15.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.