**Oral-B®**





TRIZONE powered by BRAUN

500

550

Type 3757

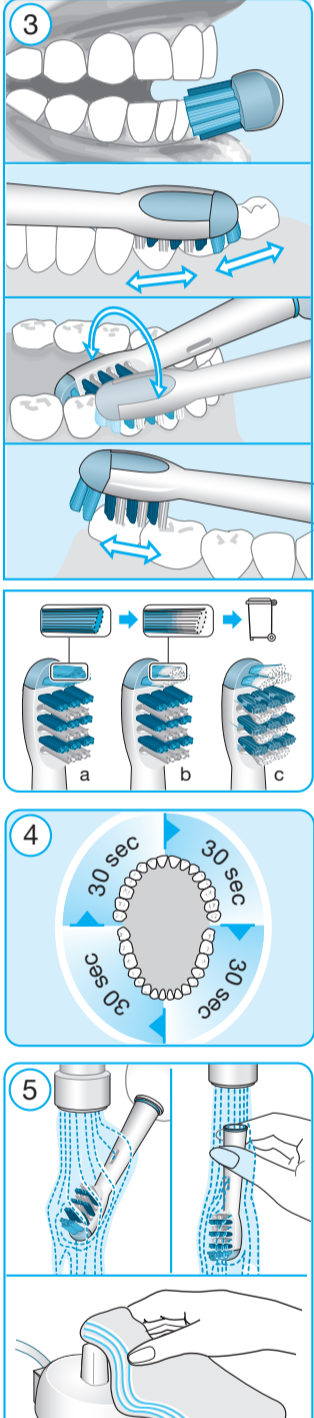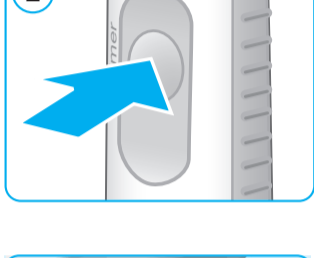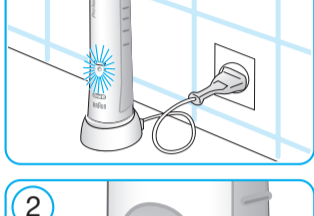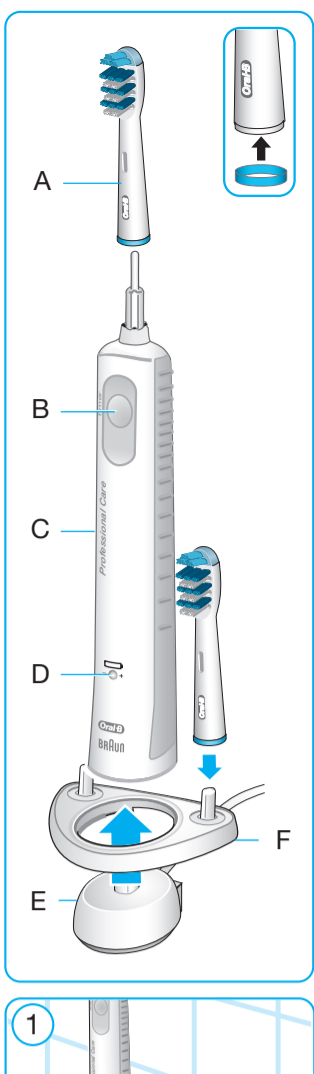- (DE) (AT) 00 800 27 28 64 63**
- (CH) 08 44 - 88 40 10**
- (GB) 0 800 783 7010**
- (IE) 1 800 509 448**
- (FR) 0 800 944 803**  
(appel gratuit depuis un poste fixe)
- (BE) 0 800 14 592**
- (LU) 800 21172**
- (ES) 901 11 61 84**
- (PT) 808 20 00 33**
- (IT) 800 124 600**
- (NL) 0800-4 45 53 88**
- (DK) 70 15 00 13**
- (NO) 22 63 00 93**
- (SE) 020 - 21 33 21**
- (FI) 020 377 877**
- (GR) 800 801 3457**

www.oralb.com  
www.braun.com  
www.service.braun.com

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## English

Before operating this toothbrush,

please read these instructions and

save this manual for future reference.

#### IMPORTANT

• Periodically check the cord for

damage. If cord is damaged, take

the charging unit to an Oral-B Service

Centre. A damaged or non-functioning

unit should no longer be used.

4

- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are supervised by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children. Children should be supervised to ensure that they do not play with the appliance.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service Centre. Do not disassemble the product except when disposing of the battery.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Do not insert any object into any opening of the appliance/charger.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

### Specification

For voltage specifications refer to the bottom of the charging unit.

### Description

- A Brush head
- B On/off button
- C Handle
- D Charge indicator
- E Basic charger
- F Brush head holder (depending on model)

### Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the basic charger (E) or the assembled charging unit (E/F) into an electrical outlet and place the toothbrush handle (C) on the charging unit.
- The green charge indicator light (D) flashes while the handle is being charged. Once it is fully charged, the green charge indicator light shines continuously for 5 seconds and then it turns off. A full charge enables up to seven days of regular brushing (twice a day, 2 minutes) (picture 1).
- If the rechargeable battery is low the red low charge indicator light (D) is flashing for some seconds when turning your toothbrush off.
- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit to maintain it at full power. Overcharging is impossible. However, for environmental reasons, Oral-B recommends to unplug the charger until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Please note: If the charge indicator light is not flashing immediately, continue charging. It will appear in 10–15 minutes.

### Using your toothbrush

#### Brushing technique

This brush is designed to work with all common brushing techniques. You can continue to use your current brushing technique, e.g. you can brush just like you would do with a manual toothbrush. You also may consult your dentist or dental hygienist about the right technique for you.

We recommend the following:  
Wet brush head and apply toothpaste (before switching the appliance on); you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on. You start and stop brushing by pressing the on/off button of the handle (picture 2).

Place the toothbrush bristles against the teeth at a slight angle towards the gumline. Apply light pressure and start brushing in back and forth motions. Brush first the outsides, then the insides. Finally, brush the chewing surfaces with the same back and forth motion throughout the mouth (picture 3). Brush for 2 minutes. Please ensure that you brush all four quadrants of your mouth equally.

#### Professional Timer

A short stuttering sound at 30-second intervals reminds you to brush equally all four quadrants of your mouth (picture 4). A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time. The elapsed brushing time is memorised, even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds the timer resets.

#### Pressure Sensor

Your toothbrush has a pressure control feature installed. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition you will hear a different sound while brushing. In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

### Brush heads

Oral-B TriZone brush head features triple zone cleaning action for outstanding plaque removal, even in between teeth.

Oral-B also offers you a variety of different Oral-B brush heads that fit your Oral-B toothbrush handle and can be used for tooth-by-tooth precision cleaning.

#### Oral-B FlossAction brush head

features micropulse bristles that allow for an outstanding plaque removal in the interdental areas.

#### Oral-B Precision Clean brush head

cups each tooth with a curved bristle trim and reaches deep between teeth at the same time.

#### Oral-B Sensitive brush head

features an extra soft bristle structure that is gentle on teeth and gums.

#### Oral-B 3DWhite brush head

features a special polishing cup to naturally whiten your teeth. Please note that children under 12 years old should not use the Oral-B 3DWhite brush head.

Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway down approximately in a 3-month period. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums (picture 3).

We do not recommend using the Oral-B FlossAction or Oral-B 3DWhite brush head with braces. You may use the Oral-B Ortho brush head, specifically designed to clean around braces and wires.

### **Warranty replacement brush heads**

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.

- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.

- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards.

Oral-B does not sell replacement brush heads or handle parts under any other brand name.

### **Cleaning Recommendations**


After brushing, rinse your brush head for several seconds under running water with the handle switched on.

Switch off the handle and remove the brush head. Rinse both parts separately under running water then wipe them dry before reassembling the toothbrush.

From time to time also clean the charging unit with a damp cloth (picture 5). The brush head holder (F) is dishwasher-safe.

Subject to change without notice.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life.

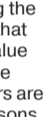Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

### **Guarantee**

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

Please also read our additional warranty statements: **Warranty replacement brush heads** on page 6.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.

### **30-Days Money Back Guarantee**

Try Oral-B TriZone starting from the day of purchase. If you are not satisfied with this product, please return your handle, refills and charging unit in their original box with your original cash register receipt to the address below within 30 days of purchase.

Great Britain/Ireland:

«Oral-B 30-Day Challenge»

Braun Cashback (BC 024)

PO Box 504

Leicester LE94 0AE

This does not affect your statutory rights. Terms and conditions apply, for further details, please contact our customer service.

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