**Oral-B®**

**Stages**  
**POWER**






powered  
by

**BRAUN**

Charger Type 3757

Handle Type 3709
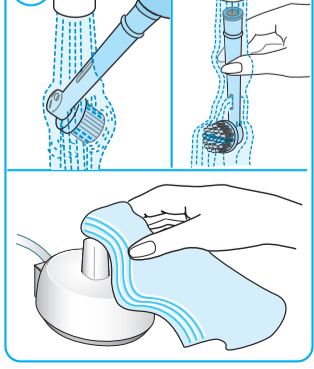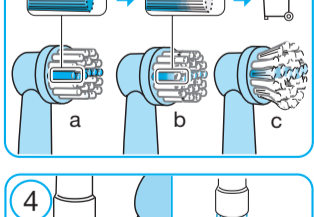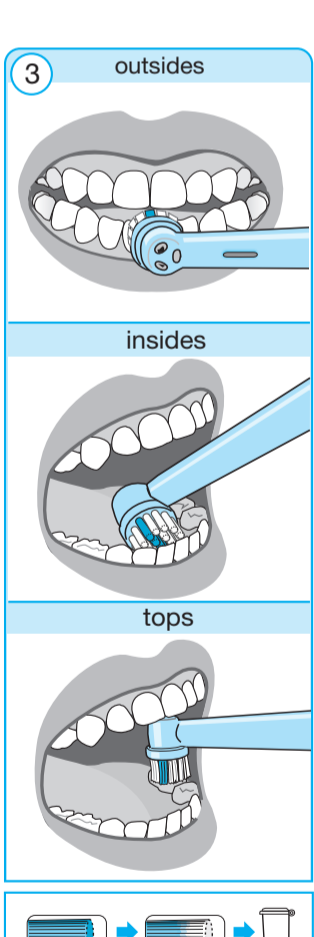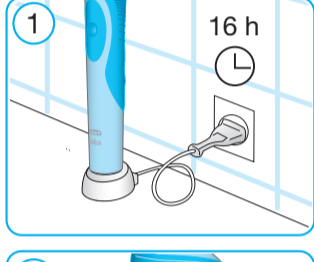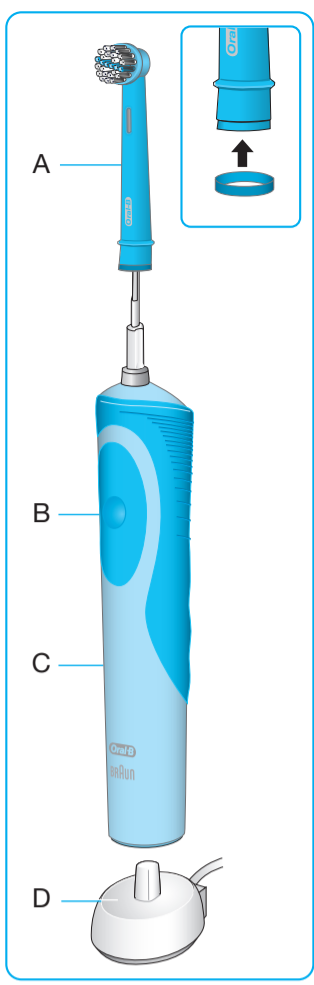
## Deutsch

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### English

*Dear parents,  
Please read the usage instructions carefully before coaching your child how to use this toothbrush. We recommend that parents supervise their children when they start using the appliance to help them to become familiar with it.*

## IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- This product is not intended for children 3 years and under. Children from 3 to 14 years, persons with reduced physical, sensory or mental capabilities, or persons with lack of experience or knowledge of the product can use the toothbrush, but should have supervision and/or instruction concerning use of the appliance in a safe way, and understand the hazards involved.
- Cleaning and maintenance should not be performed by children.
- Children should not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

### WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water, liquid or store it where it can fall or be pulled into a tub or sink. Do not reach for it if it has fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- For tooth brushing only. Do not chew or bite on brush. Do not use as a toy or without adult supervision.
- Replace brush (head) every 3 months or sooner if brush head becomes worn.
- This toothbrush is a personal care device and not intended for use on multiple patients in a dental practice or institution.

### Description

- A Brush head
- B On/off switch
- C Handle
- D Charging unit

### Specifications

For voltage specifications please refer to the bottom of the charging unit.

### Connecting and charging

This toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (D) into an electrical outlet and place the handle (C) on the charging unit ①.
- A full charge takes about 16 hours and enables up to 5 days of regular brushing (twice a day, 2 minutes).
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

### Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

### Using the toothbrush

- Wet brush head and apply kids' toothpaste. Guide the brush head to the teeth; then turn on the handle by pressing the on/off switch (B) ②.
- Start with the lower teeth: Slowly guide the brush head from tooth to tooth; first the outsides, then the insides, finally the chewing surfaces ③.

- Do not press too hard, simply let the brush do all the work.
- After two minutes of operation, an incorporated timer signals by a short stutter the end of the professionally recommended brushing time.
- Turn off the handle while the brush head is still in the mouth.

The brush head is provided with a row of blue INDICATOR® bristles ③a helping you to monitor brush head replacement need. With proper brushing twice a day for two minutes, the colour will fade halfway in a 3-month period ③b. If bristles begin to splay before the colour starts to recede, this indicates that your child exerts too much pressure on teeth and gums ③c.

#### Replacement brush heads: Oral-B Stages Power (EB10)

Available at your retailer or at an Oral-B Braun Service Centre.

In the first days of using any electric toothbrush, gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

### Cleaning

After use, rinse the brush head under running water with the handle switched on. Then switch off the handle and remove the brush head. Rinse both parts separately under running water and wipe them dry. From time to time, also clean the charging unit with a damp cloth ④.

Subject to change without notice.

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

### Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Customer Service Centre.

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