

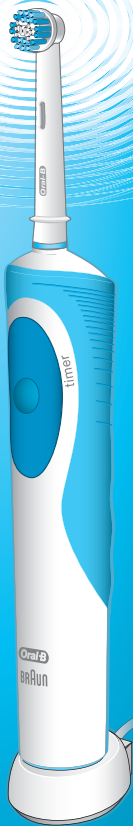
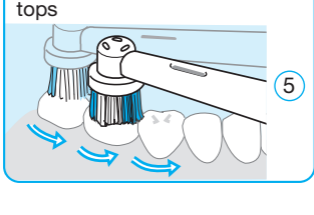
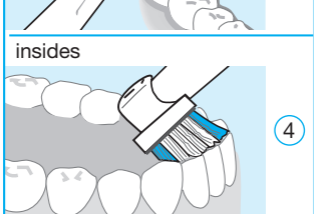
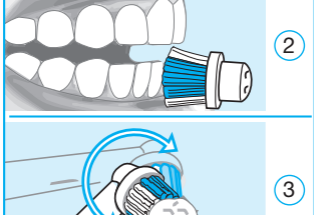
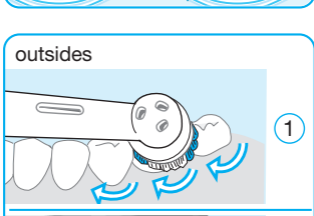
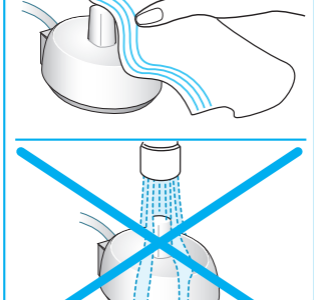
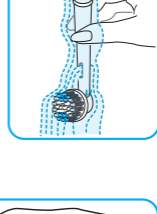
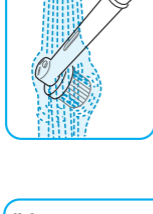
Oral-B®

VITALITY powered by BRAUN

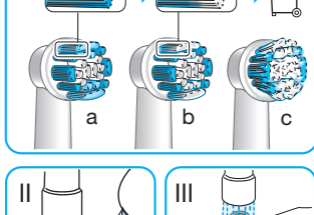Type 3757

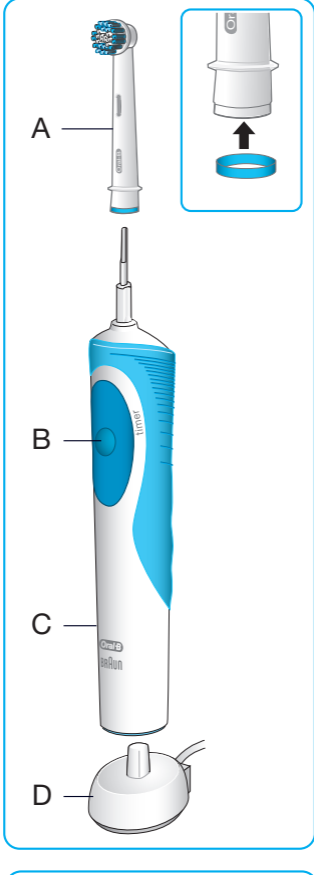- outsides**
- insides**
- tops**
- 1**
- 2**
- 3**
- 4**
- 5**
- DE AT 00 800 27 28 64 63**
- CH 08 44 - 88 40 10**
- GB 0 800 783 7010**
- IE 1 800 509 448**
- FR 0 800 944 803**
(appel gratuit depuis un poste fixe)
- BE 0 800 14 592**
- LU 800 21172**
- ES 901 11 61 84**
- PT 808 20 00 33**
- IT 800 124 600**
- NL 0800-4 45 53 88**
- DK 70 15 00 13**
- NO 22 63 00 93**
- SE 020 - 21 33 21**
- FI 020 377 877**
- GR 800 801 3457**

Internet
www.oralb.com
www.braun.com
www.service.braun.com

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English

Your Oral-B Vitality toothbrush has been developed with dental professionals to effectively remove plaque from all areas of your teeth.

Please read the use instructions carefully before first use.

IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- This product is not intended for use by children under age three. Children from 3 to 14 years, persons with reduced physical, sensory or mental capabilities, or persons with lack of experience or knowledge of the product can use the toothbrush, but should have

supervision and/or instruction concerning use of the appliance in a safe way, and understand the hazards involved.

- Cleaning and maintenance should not be performed by children.
- Children should not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water, liquid or store it where it can fall or be pulled into a tub or sink. Do not reach for it if it has fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

A Brush head (depending on model/ interchangeable)

B On/off switch

C Handle

D Charging unit

Specifications

For voltage specifications please refer to the bottom of the charging unit.

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (D) into an electrical outlet. Place the handle (C) on the charging unit.
- A full charge takes 16 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible. However for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

Using the toothbrush

The toothbrush can be used with any leading brand of toothpaste.

- Wet brush head and apply toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on (B).
- Slowly guide the brush head from tooth to tooth, spending a few seconds on each tooth surface.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces (1–5). Brush at least 2 minutes to thoroughly remove plaque. Do not press too hard or scrub.
- Turn off the handle by pressing the on/off switch.


Timer (depending on model)

An incorporated timer memorises the elapsed brushing time. A short stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

Individual brush heads and accessories

Oral-B offers you a variety of brush heads which are available in stores that sell Oral-B products. Following brush heads are interchangeable and can be used on your Oral-B Vitality handle.

Oral-B FlossAction brush head

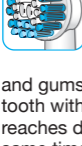
The Oral-B FlossAction brush head is recommended for everyday deep cleaning. It features micropulse bristles that get deep in hard-to-reach areas between your teeth, gently removing hidden plaque.

Oral-B Precision Clean brush head


The Oral-B Precision Clean brush head provides effective and gentle cleaning action to help protect both teeth and gums. The brush head cups each tooth with a curved bristle trim and reaches deep between teeth at the same time.

Oral-B Sensitive brush head

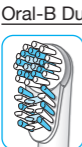
The Oral-B Sensitive brush head provides thorough cleaning for sensitive teeth and gums. It has an extra soft bristle texture that is gentle on teeth and gums.

Oral-B Dual Clean brush head

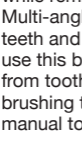The Oral-B Dual Clean brush head offers two types of cleaning motion in one. It combines the oscillating–rotating motion with an additional side-to-side motion. The oscillation motion stimulates gums while removing plaque bacteria. Multi-angled bristles clean between teeth and along the gumline. You can use this brush head either by moving from tooth to tooth or apply the same brushing technique as you do with any manual toothbrush.

Oral-B 3DWhite brush head

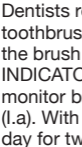For occasional or everyday polishing, please use the Oral-B 3DWhite brush head. It features a dentist-inspired polishing cup to remove surface stains from coffee, tea and tobacco for effective brightening. We do not recommend using the Oral-B 3DWhite brush head with braces.

Dentists recommend replacing any toothbrush every 3 months. Therefore the brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need (I.a). With thorough brushing, twice a day for two minutes, the blue color will fade halfway down approximately in a 3-month period (I.b). If the bristles splay before the colour recedes (I.c), you are exerting too much pressure on teeth and gums.


During the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or dental hygienist.

Cleaning

After use, rinse the brush head for several seconds under running water with the handle switched on (II). Then switch off the handle and take off the brush head. Rinse both parts separately under running water (III) and wipe them dry. From time to time, also clean the charging unit with a damp cloth (IV).

Subject to change without notice.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

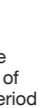Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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